VEVOR[®]

TOUGH TOOLS, HALF PRICE

Technical Support and E-Warranty Certificate www.vevor.com/support

RED LIGHT THERAPY BELT

MODEL: PLR-120S

We continue to be committed to provide you tools with competitive price.

"Save Half", "Half Price" or any other similar expressions used by us only represents an estimate of savings you might benefit from buying certain tools with us compared to the major top brands and does not necessarily mean to cover all categories of tools offered by us. You are kindly reminded to verify carefully when you are placing an order with us if you are actually saving half in comparison with the top major brands.

VEVOR®

TOUGH TOOLS, HALF PRICE

RED LIGHT THERAPY BELT

MODEL: PLR-120S

Thank you for purchasing our light, we hope you can reap benefits of red light therapy at home.


NEED HELP? CONTACT US!

Have product questions? Need technical support? Please feel free to contact us:

Technical Support and E-Warranty Certificate
www.vevor.com/support


This is the original instruction, please read all manual instructions carefully before operating. VEVOR reserves a clear interpretation of our user manual. The appearance of the product shall be subject to the product you received. Please forgive us that we won't inform you again if there are any technology or software updates on our product.

Features

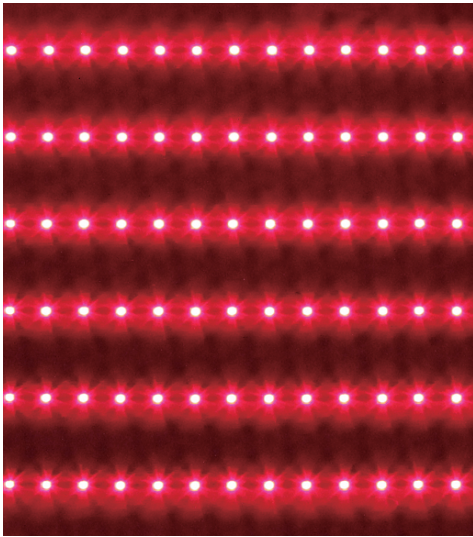
1. Perfect performance for body treatments like weight loss and pain relief.
2. Portable design so you can use it in different situations.
3. 1280 pcs leds with wavelength 660nm red light and 850nm infrared light (this light can not see by our own eyes, but they can be lighten up, you can check with you phone camera)
- 4 . We provide 1 year warranty for our pad.


NEAR INFRARED– INFRARED A (IR A) Usually use 660nm–850nm

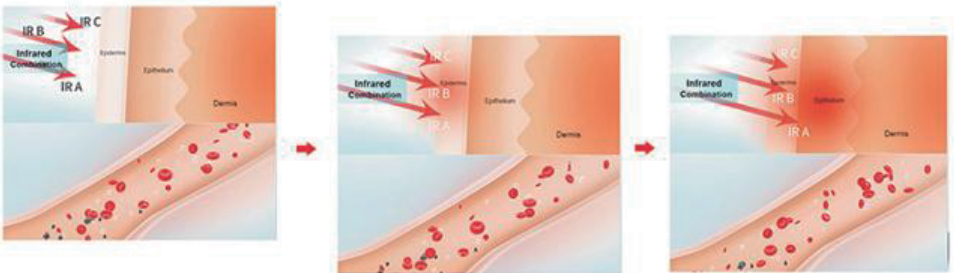Invisible near infrared energy are absorbed by muscle, subcutaneous tissue, which results in THERMAL EFFECT your body responds to this radiation heat with increased blood flow, allowing the parts of the body to receive the oxygen and nutrients they need in order to function more effectively, regeneration is stimulated.

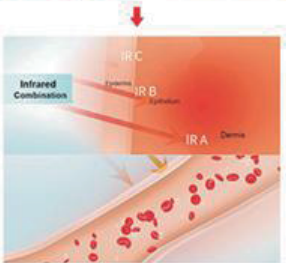Measured value

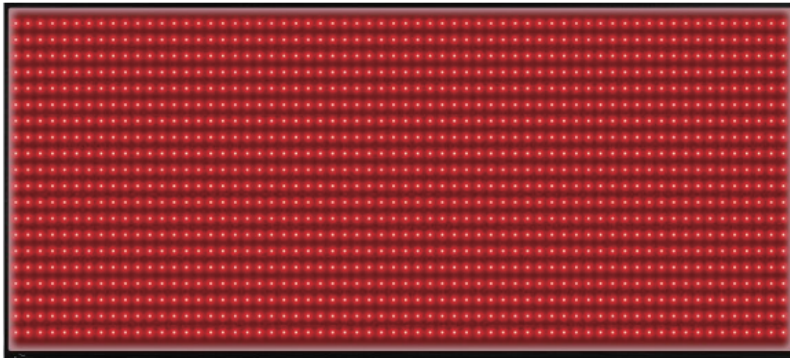Anti perspiration protective film, Sweat resistant,removable, and easy to clean.

This package includes

1*Red light Therapy Pad

1*Power Cord

1*Adaptor

1*Connecting wire

1*User Manual

1*Goggles

Specifications

Model :	Red Light therapy Pad	Wave lenath	1*660nm 2*850nm
Size :	180*80cm	Function :	Body slim/ Pain Relief
Color :	Black	Custom :	Logo/Shape/Wavelength

Power:	140W	Power Bank:	Not
--------	------	-------------	-----

Usage Frequency

Use the device daily 1-2 times and last for at least 4 weeks, and each time suggests 20 minutes on one treatment area. Then move the device to the next treatment area.

Store

Unplug the device and store it in a cool dry area.

Benefits backed by science:

Research shows , that the wavelengths in the spectrum of 580-900nm show optimal biotical response

Promote lymphocytes to metabolize toxins

Hair Re-growth

Increase collagen

Decreases Muscle and Joints pain

Enhanced blood circulation

Weight loss

Speed wound healing

Aging spots and wrinkle removal

Diminish inflammation

Red Light Energy
P1-P5

Timer
10-90 minutes options

Power
Power On, 10Hz,40Hz,Power Off

Thank you for choosing our red light therapy panel. We really appreciate your support and hope you'll enjoy this device.

SAFETY CAUTION:

1. Do not make direct view or observation to the laser beams via any optical instruments.
2. Keep out of reach of children to prevent injuries to their eyes.
3. It's forbidden for pacemaker users since the device may interfere with the operation of pacemaker.
4. During the treatment period, the diet structure should be improved and proper exercises are necessary to enhance the metabolism.
5. When using the device, please keep light emitting ports clean. Otherwise, the light losses will increase and then affect effects.
6. The symptoms like thirst, dry eyes, dizziness or rise of pulse rate after use are normal phenomena.

CORRECT DISPOSAL

This product is subject to the provision of European Directive 2012/19/EU. The symbol showing a wheeled bin crossed through indicates that the product requires separate refuse collection in the European Union. This applies to the product and all accessories marked with this symbol. Products marked as such may not be discarded with normal domestic waste, but must be taken to a collection point for recycling electrical and electronic devices.

FCC Information:

CAUTION: Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment!

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- 1) This product may cause harmful interference.

2) This product must accept any interference received, including interference that may cause undesired operation.

WARNING: Changes or modifications to this product not expressly approved by the party responsible for compliance could void the user's authority to operate the product.

Note: This product has been tested and found to comply with the limits for a Class B digital device pursuant to Part 15 of the FCC Rules, These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This product generates, uses and can radiate radio frequency energy, and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this product does cause harmful interference to radio or television reception, which can be determined by turning the product off and on, the user is encouraged to try to correct the interference by one or more of the following measures.

- Reorient or relocate the receiving antenna.
- Increase the distance between the product and receiver.
- Connect the product to an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for assistance.

VEVOR[®]

TOUGH TOOLS, HALF PRICE

Technical Support and E-Warranty Certificate

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


NEAR INFRARED- INFRARED A (IR A) Usually use 660nm-850nm

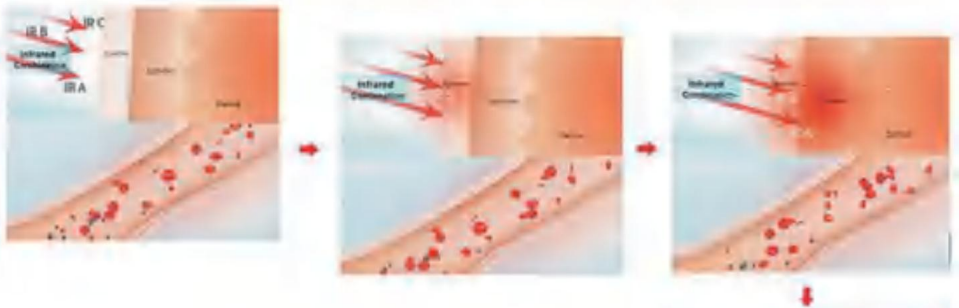invisible near infrared energy are absorbed by muscle, subcutaneous tissue, which results in THERMAL EFFECT your body responds to this radiation heat with increased blood flow, allowing the parts of the body to receive the oxygen and nutrients they need in order to function more effectively, regeneration is stimulated.

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
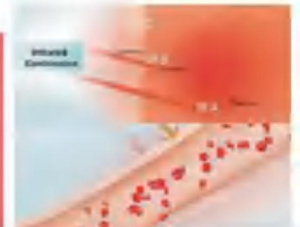
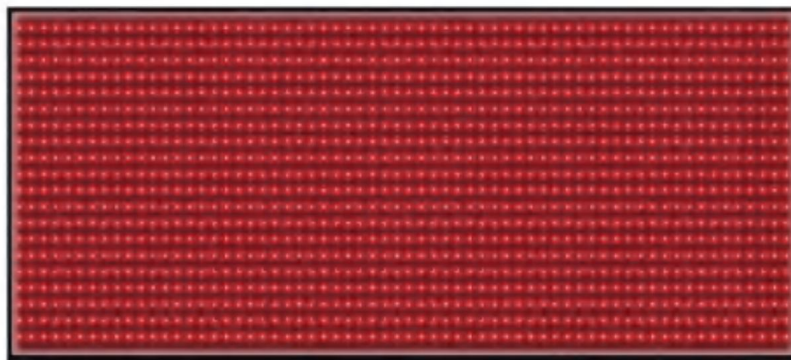
NEAR INFRARED – INFRARED A (IR A) Usually use 660nm–850nm

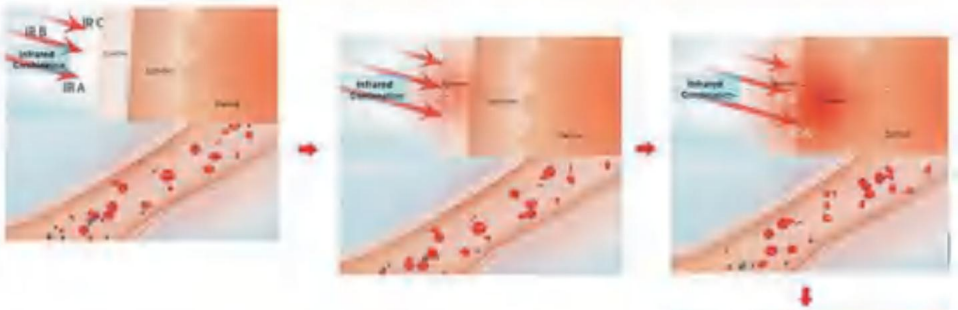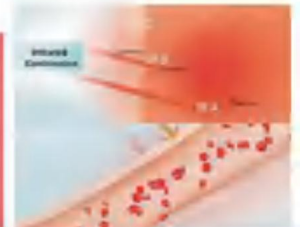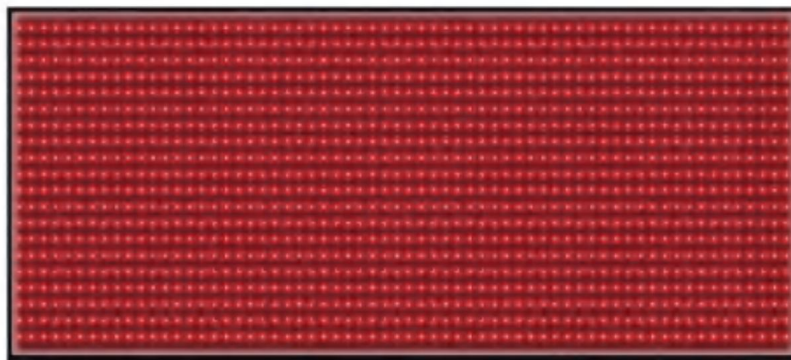invisible near infrared energy are absorbed by muscle, subcutaneous tissue, which results in **THERMAL EFFECT** your body responds to this radiation heat with increased blood flow, allowing the parts of the body to receive the oxygen and nutrients they need in order to function more effectively, regeneration is stimulated.

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