# **VEVOR<sup>®</sup>**

**TOUGH TOOLS, HALF PRICE**

Technical Support and E-Warranty Certificate [www.vevor.com/support](http://www.vevor.com/support)

## **RED LIGHT THERAPY BELT**

**MODEL: BRL-1140**

We continue to be committed to provide you tools with competitive price.

"Save Half", "Half Price" or any other similar expressions used by us only represents an estimate of savings you might benefit from buying certain tools with us compared to the major top brands and does not necessarily mean to cover all categories of tools offered by us. You are kindly reminded to verify carefully when you are placing an order with us if you are actually saving half in comparison with the top major brands.

# VEVOR®

TOUGH TOOLS, HALF PRICE

## RED LIGHT THERAPY BELT

MODEL: BRL-1140

### NEED HELP? CONTACT US!


Have product questions? Need technical support? Please feel free to contact us:

**Technical Support and E-Warranty Certificate**  
[www.vevor.com/support](http://www.vevor.com/support)

This is the original instruction, please read all manual instructions carefully before operating. VEVOR reserves a clear interpretation of our user manual. The appearance of the product shall be subject to the product you received. Please forgive us that we won't inform you again if there are any technology or software updates on our product.

# Instructional Manual


## NEAR INFRARED – INFRARED A (IR A) Usually use 660nm–850nm

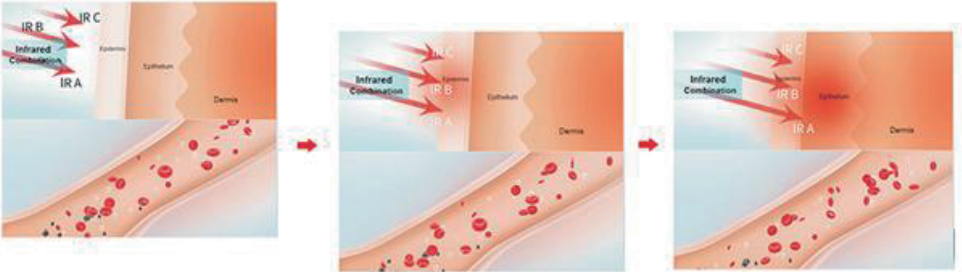Invisible near infrared energy are absorbed by muscle, subcutaneous tissue, which results in THERMAL EFFECT your body responds to this radiation heat with increased blood flow, allowing the parts of the body to receive the oxygen and nutrients they need in order to function more effectively, regeneration is stimulated.


## Measured value


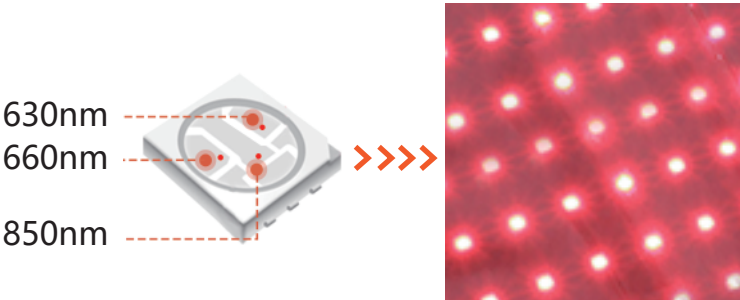
## Features:

- 1: Perfect performance for body treatments like weight loss and pain relief.
- 2: Portable design so you can use it in different situations.
3. 264 pcs leds with wavelength 630nm 660nm red light and 850nm infrared light ( this light can not see by our own eyes, but they can be lighten up, you can check with you phone camera).
4. Please press the on&off button for 3 seconds when using the switch.
5. We provide 1year warranty for our belt.

## weight loss and pain relief
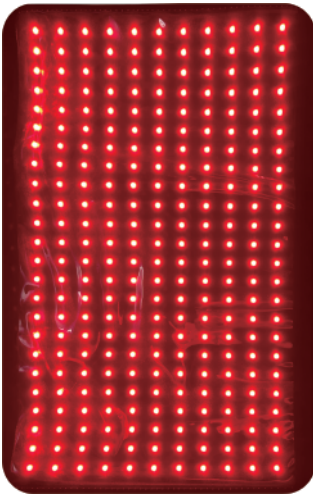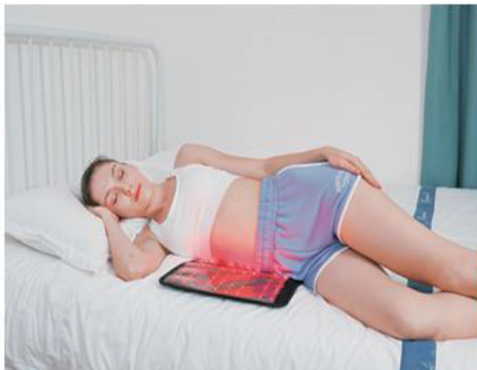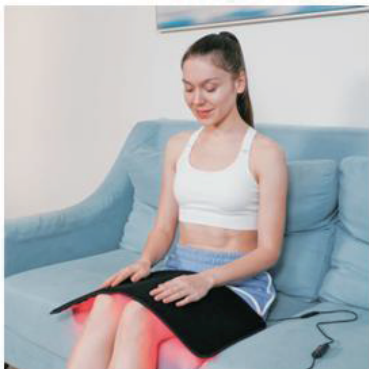
## Specifications:

Rating power	132W
Actual Power	26W
Product Size	550*350*6 MM
LED QTY	264*0.5W
Net Weight	1.25KG
Wavelength	1*630nm;1*660nm;1*850nm
EMF	OUT
Irradiance	Max: 40mW/cm <sup>2</sup>
Timer	M5/10/15/20/25/30min
Material	Fabric
Lifespan	Over 50,000 Hours

### This package includes:

- 1\* Red light therapy pad
- 1\* Power cord
- 1\* Adaptor
- 1\*Goggles

- 1\* Controller
- 2\* Stretchy straps
- 1\* User manual

## SAFETY CAUTION:

1. Do not make direct view or observation to the laser beams via any optical instruments.

2. Keep out of reach of children to prevent injuries to their eyes.
3. It's forbidden for pacemaker users since the device may interfere with the operation of pacemaker.
4. During the treatment period, the diet structure should be improved and proper exercises are necessary to enhance the metabolism.
5. When using the device, please keep light emitting ports clean. Otherwise, the light losses will increase and then affect effects.
6. The symptoms like thirst, dry eyes, dizziness or rise of pulse rate after use are normal phenomena.

## CORRECT DISPOSAL

This product is subject to the provision of European Directive 2012/19/EU. The symbol showing a wheeled bin crossed through indicates that the product requires separate refuse collection in the European Union. This applies to the product and all accessories marked with this symbol. Products marked as such may not be discarded with normal domestic waste, but must be taken to a collection point for recycling electrical and electronic devices.

## FCC Information:

**CAUTION:** Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment!

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- 1) This product may cause harmful interference.
- 2) This product must accept any interference received, including interference that may cause undesired operation.

**WARNING:** Changes or modifications to this product not expressly approved by the party responsible for compliance could void the user's authority to operate the product.

**Note:** This product has been tested and found to comply with the limits for a Class B digital device pursuant to Part 15 of the FCC Rules, These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This product generates, uses and can radiate radio frequency energy, and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this product does cause harmful interference to radio or television reception, which can be determined by turning the product off and on, the user is encouraged to try to correct the interference by one or more of the following measures.

- Reorient or relocate the receiving antenna.
- Increase the distance between the product and receiver.
- Connect the product to an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for assistance.

**VEVOR<sup>®</sup>**

**TOUGH TOOLS, HALF PRICE**

**Technical Support and E-Warranty Certificate**

**[www.vevor.com/support](http://www.vevor.com/support)**

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

### NEAR INFRARED – INFRARED A (IR A) Usually use 660nm–850nm

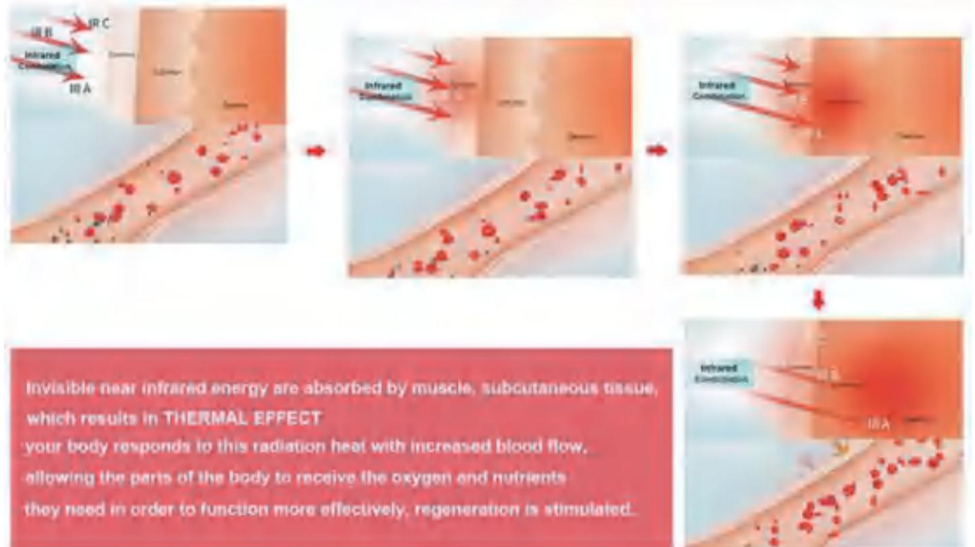Invisible near infrared energy are absorbed by muscle, subcutaneous tissue, which results in **THERMAL EFFECT** your body responds to this radiation heat with increased blood flow, allowing the parts of the body to receive the oxygen and nutrients they need in order to function more effectively, regeneration is stimulated.

The diagram illustrates the process of near-infrared (IR A) absorption. It shows three stages: 1. IR B and IR C rays are shown being reflected or passing through, while IR A rays are absorbed by the skin. 2. The absorption of IR A energy causes the skin and underlying tissues to heat up. 3. This heating leads to a thermal effect, which increases blood flow in the capillaries, bringing oxygen and nutrients to the area. A text box explains that this process stimulates regeneration.

## Messwert

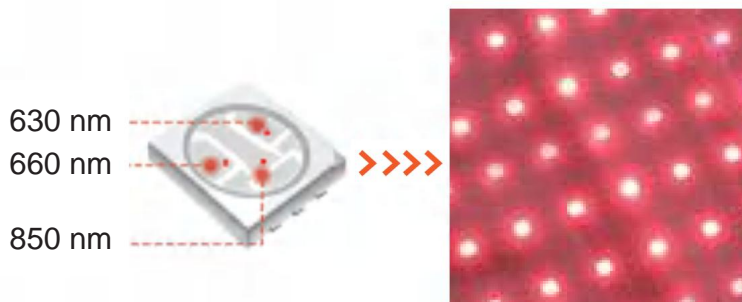630 nm  
660 nm  
850 nm

The diagram shows a device with three red dots indicating light emission at 630 nm, 660 nm, and 850 nm. Four orange arrows point from the device to a square image of red light particles, representing the emitted near-infrared radiation.

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
The diagram illustrates the process of infrared energy absorption. It shows a cross-section of skin and muscle layers. Red arrows labeled 'Infrared Contribution' point towards the tissue. Below the skin, red dots represent blood cells in a vessel, with arrows indicating increased blood flow. A text box explains that this leads to a thermal effect and stimulates regeneration.

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
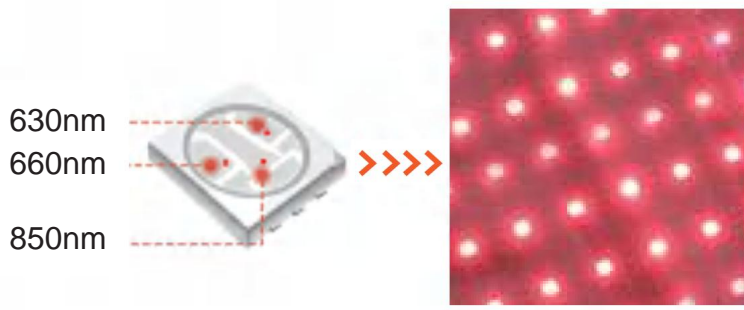
**NEAR INFRARED – INFRARED A (IR A)**  
Usually use 660nm–850nm

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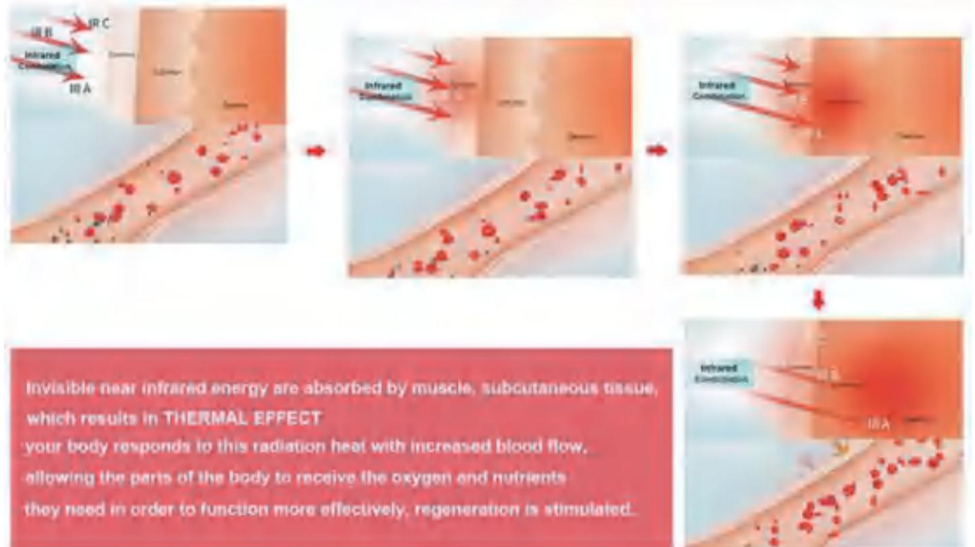The diagram illustrates the process of near-infrared (IR A) absorption. It shows three stages: 1) IR B and IR C rays are blocked by the skin, while IR A rays penetrate. 2) IR A rays reach the muscle and subcutaneous tissue, where they are absorbed. 3) This absorption causes a thermal effect, leading to increased blood flow and the delivery of oxygen and nutrients to the tissue.

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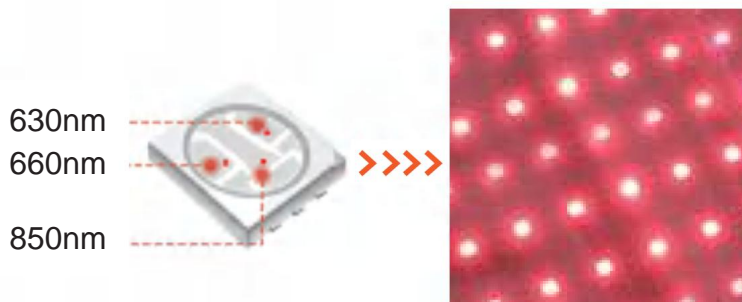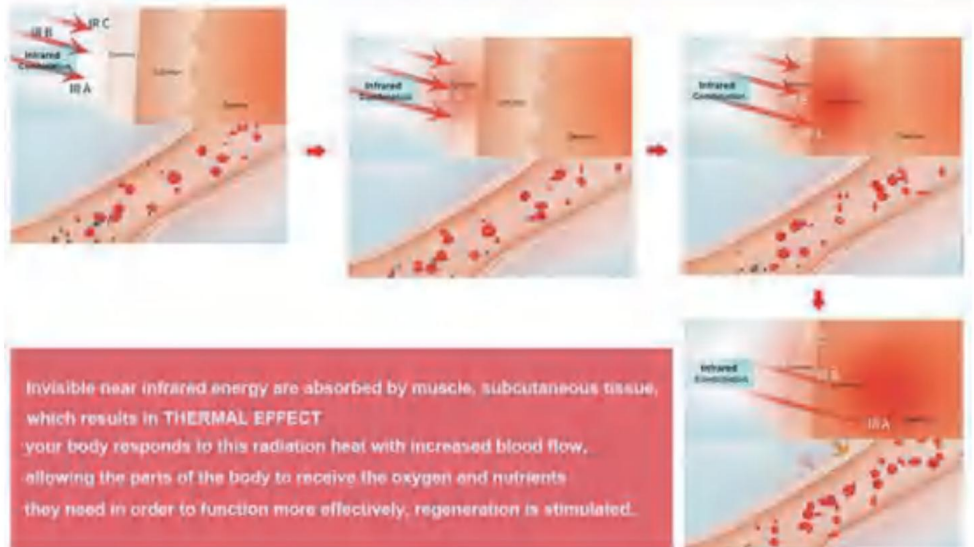
The diagrams show a cross-section of skin with layers labeled: Epidermis, Dermis, and Muscle. Infrared radiation (IR A) is shown entering from the left. The first diagram shows the radiation passing through the epidermis and dermis into the muscle. The second diagram shows the radiation being absorbed by the muscle, causing it to heat up. The third diagram shows the heat causing blood vessels in the muscle to dilate, increasing blood flow. The fourth diagram shows the increased blood flow bringing oxygen and nutrients to the muscle, stimulating regeneration.

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