# **VEVOR<sup>®</sup>**

## **TOUGH TOOLS, HALF PRICE**

Technical Support and E-Warranty Certificate [www.vevor.com/support](http://www.vevor.com/support)

### **LED FACE THERAPY**

**MODEL: FM60;FM60X**

We continue to be committed to provide you tools with competitive price.

"Save Half", "Half Price" or any other similar expressions used by us only represents an estimate of savings you might benefit from buying certain tools with us compared to the major top brands and does not necessarily mean to cover all categories of tools offered by us. You are kindly reminded to verify carefully when you are placing an order with us if you are actually saving half in comparison with the top major brand

# VEVOR®

TOUGH TOOLS, HALF PRICE

## LED FACE THERAPY

Model: FM60

FM60X

### NEED HELP? CONTACT US!

Have product questions? Need technical support? Please feel free to contact us:

**Technical Support and E-Warranty Certificate**  
**[www.vevor.com/support](http://www.vevor.com/support)**

This is the original instruction, please read all manual instructions carefully before operating. VEVOR reserves a clear interpretation of our user manual. The appearance of the product shall be subject to the product you received. Please forgive us that we won't inform you again if there are any technology or software updates on our product.

# IMPORTANT USE PRECAUTIONS

## WARNING:

1. This product is not a toy. Do not allow children to play with or approach this item.
2. Use as intended only. Do not stand on the product.
3. TO AVOID ELECTRICAL SHOCK, DO NOT put liquid of any kind into the Main Unit Housing containing the electrical components.
4. When installing and using, check whether the nails installed on the wall are firm
5. For your safety, please do not disassemble or modify the product

## SAFETY INSTRUCTIONS, WARNINGS, AND DISCLAIMER

### SUPERVISION

1. **DO NOT** use this device on infants.
2. **DO NOT** let sleeping, disabled, or unconscious persons use the device.
3. **DO NOT** exceed 20 minutes of use in one therapy session.
4. **DO NOT** use on open, fresh wounds.
5. **DO NOT** use in combination with liniments, salves, ointments, or balms that contain heat-producing ingredients. A skin burn could be the result of doing so.
6. **DO NOT** modify or attempt to repair this equipment. There are no parts that are serviceable by the user.

**This device is an electrical device. To avoid electric shock and other electricity-related dangers, adhere to the following instructions:**

1. **DO NOT** wash electrical parts with water or other liquids.
2. **DO NOT** apply or touch the device to, or with, wet or damp skin.
3. **DO NOT** pull, carry, or lift the device by its cord. If the cord is damaged, don't use the device.
4. **DO NOT** use on or near heated surfaces.
5. **DO NOT** operate this device in areas where it could be exposed to flammable or combustible products or vapours. Explosion of fire may occur.
6. **DO NOT** store your device in hot temperatures exceeding 140 degrees Fahrenheit.
7. This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction

concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.

**The product is constantly being upgraded, and there will be no notification if there is any change**

**SAVE THESE INSTRUCTIONS**

## **FCC Information**

### **CAUTION:**

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment!

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- 1) This product may cause harmful interference.
- 2) This product must accept any interference received, including interference that may cause undesired operation.

### **WARNING:**

Changes or modifications to this product not expressly approved by the party responsible for compliance could void the user's authority to operate the product.

### **Note:**

This product has been tested and found to comply with the limits for a Class B digital device pursuant to Part 15 of the FCC Rules, These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This product generates, uses and can radiate radio frequency energy, and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this product does cause harmful interference to radio or television reception, which can be determined by turning the product off and on, the user is encouraged to try to correct the interference by one or more of the following measures.

- Reorient or relocate the receiving antenna.
- Increase the distance between the product and receiver.
- Connect the product to an outlet on a circuit different from that to which the receiver is connected.

- Consult the dealer or an experienced radio/TV technician for assistance.

## **THE BENEFITS OF RED LIGHT THERAPY FOR DIVERSE CONDITIONS**

A few of the scientifically proven results of red light therapy investigated thus far include:

### **1. Reduced Inflammation**

Near-infrared light reduces overall inflammation in the body by decreasing the presence of inflammatory markers. This reduction can help ease the symptoms associated with joint pain, sore muscles, autoimmune diseases, arthritis, traumatic brain injuries, and spinal cord injuries.

### **2. Reduced Pain**

In one clinical study, neuropathic pain caused by a spinal cord injury was dramatically reduced by the application of red light treatment. Scientists conducting the study referred to it as “an exciting prospect with significant clinical relevance.” Near-infrared light wavelengths reduce overall pain by easing joint stiffness and soreness, diminishing inflammation, easing muscle spasms and enhancing blood flow.

### **3. Reduced recovery times after training and injury**

For high-performance athletes (and those who train like them) red light therapy can help accelerate muscle repair following fatigue and injury. Mitochondria within cells are particularly responsive to red light therapy, and muscle cells are exceptionally rich in mitochondria. Red light therapy may additionally stimulate stem cells, further assisting in muscle recovery.

#### **Healthier, plumper-looking skin due to increased collagen production.**

Red light therapy can dramatically transform the skin.

Red light wavelengths in particular target the mitochondrial chromophores within skin cells, generating the production of collagen proteins. Collagen stimulation yields more holistic and enduring benefits than simply resurfacing the outer layers of the skin. Stem cells may also be activated, increasing tissue repair. The result is accelerated healing and wound repair, improved appearance in hypertrophic scars, a reduction in fine lines and wrinkles, and improved skin texture.

Participants in one study experienced a clinically assessed improvement in skin complexion and collagen density. Improved circulation (improved blood flow). In many studies, red light therapy has been clinically proven to increase the diameter of blood vessels and to improve circulation. Additionally, red light therapy protects red blood cells against oxidative stress and limits platelet loss during surgical procedures.

#### **4. What types of light are there and what do they do?**

Antiaging, rejuvenation, and perfecting your complexion: fine lines, loss of elasticity, dehydration, large pores, dull, uneven tone or sun-damaged skin, use combined red and NIR light on your face, neck, chest, and hands.

Heals & relieves:

- ✓ Psoriasis
- ✓ Muscle & joint pain
- ✓ Thins & delicate skin, wounds, burns & red scars

To regulate sebum and reduce and control mild to moderate acne, add blue light to the combo.

#### **Red & NIR**

Antiaging, rejuvenation, and perfecting your complexion: fine lines, loss of elasticity, dehydration, large pores, dull, uneven tone or sun-damaged skin, use combined red and NIR light on your face, neck, chest, and hands.

Rejuvenates & perfects:

- ✓ Brightens & hydrates
- ✓ Firms, tones & depuffs
- ✓ Smooths bumps & wrinkles
- ✓ Clarifies redness & refines pores
- ✓ Fades fine lines, spider veins, dark spots & scars

And use red and NIR again for psoriasis and wound healing such as sunburn, open wounds and red scars, thin and vulnerable skin, and muscular or joint pain relief on your face, body, and scalp.

Heals & relieves:

- ✓ Psoriasis
- ✓ Muscle & joint pain
- ✓ Thins & delicate skin, wounds, burns & red scars

#### **Blue, Red & NIR**

To regulate sebum and reduce and control mild to moderate acne, add blue light to the combo.

Mild to moderate acne:

- ✓ Kills acne causing bacteria
- ✓ REgulates sebum/ oily T-zone
- ✓ Controls & heals blemishes, spots & comedones

#### **NIR**

Use NIR for dry, tight, flaky, and rough skin, redness such as diffused redness and visible capillaries, or compromised, thin and vulnerable skin.

Sensitive skin:

- ✓ Diffused redness
- ✓ Visible capillaries
- ✓ Dry, tight, flaky, or rough skin
- ✓ Compromised, thin or vulnerable skin

#### **Red**

Use red for hyperpigmentation (dark spots) on your face caused by medicines, hormones, and post inflammatory hyperpigmentation (PIH)

- Hyperpigmentation & dark spots caused by:
- ✓ Medicines
  - ✓ Hormones
  - ✓ Post inflammatory hyperpigmentation (PIH)

## Product Introduction

### FM60

#### Technical Specifications

**LED Spectra:** 630nm&660nm&850nm

**Color:** Black

**Number of LEDs:** 222LEDs (3per LED bulb)  
74LED bulbs

**Face Mask Size:** 300mm(w)x195mm(h)

**Program time:** 10minutes

**Electrical output:** 5V/1A/2500mAh

**Operating Temperature:** 0°C to 40°C

**Storage and Transport Temperature:** -5°C to 55°C


#### Model

- 1.High-end
- 2.Mid-range
- 3.Low-end

#### Device Description

1. Silicon flexible face mask
2. Controller
3. Power supply and country specific adaptors
4. USB C to USB A connector
5. Head straps
6. Carry bag
7. Charger
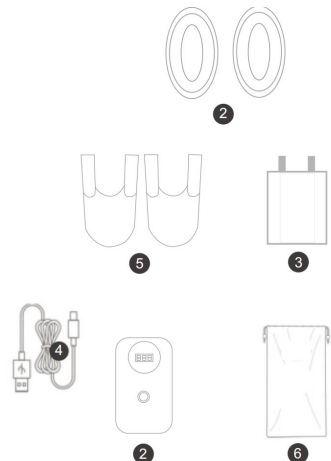
# FM60X

## Technical Specifications

**LED Spectra:** 630nm&660nm&850nm Or customized wavelengths

**Color:** Black

**Material:** Silicone rubber

**Number of LEDs:** 222LEDs(3per LED bulb)  
74LED bulbs


**Face Mask Size:** 300mm(w)x195mm(h)

**Program time:** 10minutes

**Electrical output:** 5V/1A/2500mAh


**Operating Temperature:** 0°Cto40°C

**Storage and Transport Temperature:** -5°C to 55°C
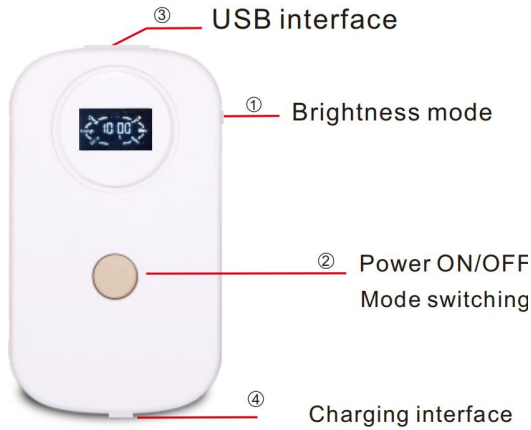
## Device Description

1. Silicon flexible face mask
2. Controller
3. Power supply and country specific adaptors
4. USB C to USB A connector
5. Head straps
6. Carry bag
7. Charger
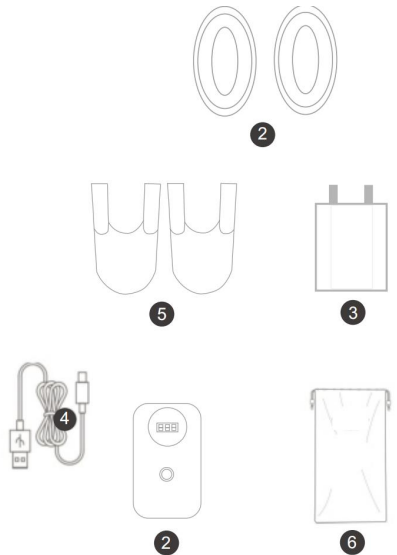
# Contraindications

Stop using the This device facial mask if you encounter any prolonged skin redness or discomfort that lasts for more than 2 hours.

If a problem occurs for more than 2 hours you might be light sensitive, or have another condition. We recommend you consult with a medical professional for advice.

1. Please do not use this product if you have lupus erythematosus, photosensitive eczema or leukemia. Otherwise, it may cause severe skin reactions.
2. Please do not use this product if you are sensitive or dizzy to light. Otherwise, it may cause severe skin reactions. Contraindications
3. It is recommended not to use this product if you have migraine.
4. Please do not use this product if you have a history of seizures caused by light

# Warnings

Safety warning:

1. Please ensure the USB port been linked to the plug adapter before plugging in the power source.
2. Do not place the device at spots of water or dripping water. Please unplug the power and stop using immediately if the device has liquid or water on.
3. If the mask light is scratched or intruded by liquid, please stop using it immediately.
4. Do not use the mask when in the sauna or in the rain.
5. Keep the mask or wires away from heat sources.
6. If the mask, controller or wire is damaged, please do not use it again.
7. Do not pull the USB cable forcefully, in case the mask cable fall off.
8. Please stop use the mask when it is overheated.
9. Do not let children use the mask without supervise.

# Care and Maintenance

Store This Device (FACE MASK) system at room temperature in dry conditions. Use the carry bag provided.

Before each use, with the device switched off, clean the facemask following the cleaning instructions below.

Avoid storing This Device LED system in dusty conditions, in direct sunlight or in damp and humid conditions.

# COEEECT DISPOSAL

This product is subject to the provision of European Directive 2012/19/EU. The symbol showing a wheeled bin crossed through indicates that the product requires separate refuse collection in the European Union. This applies to the product and all accessories marked with this symbol. Products marked as such may not be discarded with normal domestic waste, but must be taken to

a collection point for recycling electrical and electronic devices.

# BATTERY DISPOSAL

Li-Ion

To minimize hazards to health and the environment at the end of this product's life, laws dealing with Waste Electrical and Electronic Equipment (WEEE) and The Waste Battery Directive require you to

dispose of this product at a suitable collection facility where it will be sent in order to remove the batteries and for appropriate recycling. Please contact your local authorities for more details on recycling and safe disposal of these in your area.

**Manufacturer:** Shanghaimuxinmuyeyouxiangongsi

**Address:** Shuangchenglu 803nong11hao1602A-1609shi, baoshanqu, shanghai 200000 CN.

**Imported to AUS:** SIHAO PTY LTD.

1 ROKEVA STREET EASTWOOD NSW 2122 Australia

**Imported to USA:** Sanven Technology Ltd.

Suite 250, 9166 Anaheim Place, Rancho Cucamonga, CA 91730

UK	REP
----	-----

YH CONSULTING LIMITED.

C/O YH Consulting Limited Office 147,  
Centurion House, London Road, Staines-  
upon-Thames, Surrey, TW18 4AX

EC	REP
----	-----

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Mainzer Landstr.69,

60329 Frankfurt am Main.

**VEVOR<sup>®</sup>**

**TOUGH TOOLS, HALF PRICE**

**Technical Support and E-Warranty Certificate**  
**[www.vevor.com/support](http://www.vevor.com/support)**

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**Importerad till AUS:** SIHAO PTY LTD.
**Importerad till USA:** Sanven Technology Ltd.
E-CrossStu GmbH