VEVOR[®]

TOUGH TOOLS, HALF PRICE

Technical Support and E-Warranty Certificate www.vevor.com/support

LIGHT THERAPY FOR FOOT

MODEL: PLR-100 PLR-100S

We continue to be committed to provide you tools with competitive price.

"Save Half", "Half Price" or any other similar expressions used by us only represents an estimate of savings you might benefit from buying certain tools with us compared to the major top brands and does not necessarily mean to cover all categories of tools offered by us. You are kindly reminded to verify carefully when you are placing an order with us if you are actually Saving Half in comparison with the top major brands.

VEVOR®

TOUGH TOOLS, HALF PRICE

LIGHT THERAPY
FOR FOOT

MODEL: PLR-100 PLR-100S

NEED HELP? CONTACT US!

Have product questions? Need technical support? Please feel free to contact us:

Technical Support and E-Warranty Certificate
www.vevor.com/support

This is the original instruction, please read all manual instructions carefully before operating. VEVOR reserves a clear interpretation of our user manual. The appearance of the product shall be subject to the product you received. Please forgive us that we won't inform you again if there are any technology or software updates on our product.

CORRECT DISPOSAL

This product is subject to the provision of European Directive 2012/19/EU. The symbol showing a wheeled bin crossed through indicates that the product requires separate refuse collection in the European Union. This applies to the product and all accessories marked with this symbol. Products marked as such may not be discarded with normal domestic waste, but must be taken to a collection point for recycling electrical and electronic devices.

MANUAL

Before using this product please read this user instruction carefully and keep it for future reference.

Safety requirements and Notes

Dear Users,

Thanks for purchasing red light & infrared light therapy. instruction for use has been provided to ensure that you will be able to understand and safely use this device.

Important Notice

If you have any problem or need help in using, please contact our technical service center for help in time. We will give you technical support or arrange professional technical expert for service at the first time.

Use correctly can extend the life of equipment, also make full usability of equipment by farthest.

Abnormal operation, hazards to equipment or personal injury caused by operations that are absolutely prohibited as specified in this document. Disclaims any responsibility for safety, reliability or performance of this equipment by not observing the instructions! Ant

faults arising from such non-observance will invalidate the warranty!
Please carefully read the instructions in the “safety requirements”, “Note” and the special warnings”! part of the content.

1. Relax muscles and relieve muscle cramps, pain and stiffness
2. Relieve foot pain
3. Accelerate healing, relieve pain and relax muscles.
4. Reduce pain and stiffness
5. Promote blood circulation
6. Enhance recovery ability after sports or injuries.


Introduction

Effective wavelength of near infrared light therapy include 750nm to 1,400nm. The wave length of the near infrared is more beneficial in treating deep wounds and relieving muscle and joint pain. 850nm is common for near infrared light therapy. Shoes adopt 66pcs 850nm infrared led lights(invisible to naked light) and 66pcs 660nm red lights. The combination of red and

infrared light increase blood circulation by thermal effect to provide pain relief.

How does infrared light work?

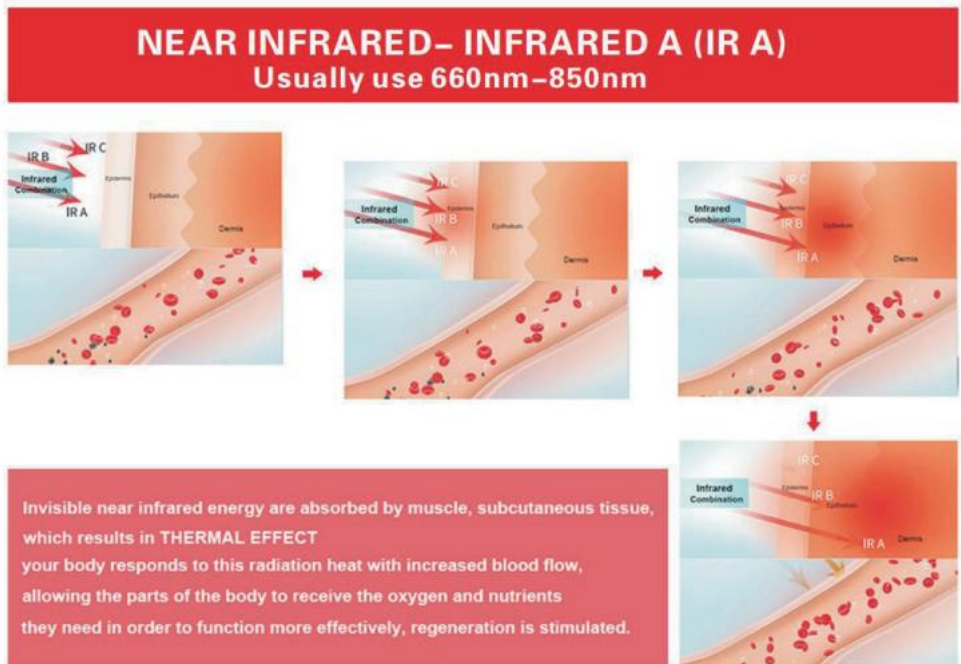The NIR is the shortest wavelength and the most penetrating infrared ray. Which can penetrate 5 to 10mm of the human body tissues.

Infrared light can penetrate the skin, directly to the muscles, subcutaneous tissue which helps to accelerate blood circulation, increase metabolism, relieve pain and increase muscle relaxation, resulting in massage effects. Infrared treats the disease at different levels because of mobilize the body's own disease-resistant ability.

► **Function**

The Amazing Light Therapy Product (infrared LED therapy Mitten) benefits:

- ◆ Relieves stiffness & spasm
- ◆ Reduces pain and inflammation
- ◆ Easily administer at home
- ◆ Improves hydration and oxygenation of the tissues
- ◆ Improve circulations
- ◆ Safe for all types of skin
- ◆ Calms the nervous system

SAFE & NATURAL PAIN TREATMENT


- ◆ Double sides treatment
- ◆ OTC class II medical device
- ◆ Relieves pain and stiffness
- ◆ Increase blood circulation
- ◆ Enhances recovery from exercise or injury
- ◆ Accelerates healing, soothes pain, and relaxes muscles

Medical -grade LED lights reach deep into joints and soft tissue Infrared & Red Light Therapy Device penetrates muscles and joints effectively to improve blood circulation and relieve pain, suitable for hand joint rheumatism nerve tendon sheath pain, etc

Targeting The Body to Protect Your feet

Running for a long time, often wearing high heels, often play football Causes manual work, and strenuous exercise leads to muscle strain.

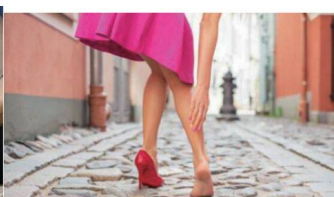Product description

LEd red light consists of infrared light shoe120pcs Lights: SMD(60pcs red lights&60pcs IR lights)

- DC cord
- AC Adapter: AC100~240V with a cord.

Instruction for use

- ◆ A1: Red&Infrared LED Lights Therapy Devices
- ◆ A2: DC Cord.
- ◆ A3:Auto Shut-Off Timer(30 Minutes)
- ◆ A4:Adapter(AC100-240V)& Power Cord.

► Cleaning

1. Disconnect the device from the power supply before cleaning.
2. Using a soft cloth with a little water to gently wipe the surface of device.

► Caution

1. Do not use alcohol or other corrosive cleaning agents to avoid corroding the surface and damaging internal parts.
2. Do not apply liquid cleaner directly to any surface of the adapter, timer, pad as it could cause damage to the electrical components.

► Storage

The device must be stored in dry area; away from high humidity

Warning

No washing

No compression

Stay dry

Handle lightly

Attention

People with heart disease or heart pacemaker can not use this light therapy system Pregnant woman, skin-ulcers and allergic people can not use it.

Features

Model:	
Material:	75% neoprene and 25%velvet
Power supply	AC 100-240Vor DC 5 V Portable power
Product size	12.5*5.3*2.3inches
wavelength	660nm+850nm
Temperature	Max 104 °F
Frequency	50/60 Hz
N.W	0.80kg
G.W	0.82kg

Warranty and disclaimer

Limited One-Year Warranty

We service a limited one-year warranty for the original purchaser against any technical defects

Warranty coverage applies to the original owner and to the original product only and is not transferable.

We warranty that the device shall be free from defects in material and workmanship for one year from the date of purchase when the device is used and maintained according to the requirements outlined in the User Manual.

If your device stops functioning or working as intended during the one-year warranty period please contact with us and one of our customer service team members will work with you to understand the issue and give you instructions on how to get your LIGHT THERAPY repaired or replaced, free of charge.

This warranty does not cover normal wear of the device and does not apply to any device that has been tampered with or damaged due to misuse, abuse, or negligent handling. Dear Customer, thanks so much for your great support on us! if there is anything problem with device, please feel free to contact us and we will be glad to help you.

Please send us the photos or videos of the product failure, and please kindly advise me your order number, phone number, contact name and shipping address.

If your order is still on one-year warranty, we will resend you the new replacement once we confirm the exact problems. Thanks for your cooperation in advance.

Specific Warranty Policy

-  1-Year defects liability period
-  60 days return and replacement warranty
-  18-month free replacement for charger, controller and connectors.
(if they are broken).
-  Return postage borne by customers

- ✓ Return requests as "I don't like it", "i don't want it", "the color doesn't look good" are not accepted.
- ✓ The less-than-three-light-broke-situation will not be covered in the warranty policy because it doesn't affect usage effect.
- ✓ keep the original outer package for successful return.

FCC Information

CAUTION: Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment!

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- 1) This product may cause harmful interference.
- 2) This product must accept any interference received, including interference that may cause undesired operation.

WARNING: Changes or modifications to this product not expressly approved by the party responsible for compliance could void the user's authority to operate the product.

Note: This product has been tested and found to comply with the limits for a Class B digital device pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This product generates, uses and can radiate radio frequency energy, and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this product does cause harmful interference to radio or television reception, which can be determined by turning the product off and on, the

user is encouraged to try to correct the interference by one or more of the following measures.

- Reorient or relocate the receiving antenna.
- Increase the distance between the product and receiver.
- Connect the product to an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for assistance

Manufacturer: Shanghaimuxinmuyeyouxiangongsi

Address: Shuangchenglu 803nong11hao1602A-1609shi, baoshanqu, shanghai 200000 CN.

Imported to AUS: SIHAO PTY LTD. 1 ROKEVA STREETEASTWOOD NSW 2122 Australia

Imported to USA: Sanven Technology Ltd. Suite 250, 9166 Anaheim Place, Rancho Cucamonga, CA 91730

UK	REP
-----------	------------

YH CONSULTING LIMITED. C/O YH Consulting Limited Office 147, Centurion House, London Road, Staines-upon-Thames, Surrey, TW18 4AX

EC	REP
-----------	------------

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VEVOR[®]

TOUGH TOOLS, HALF PRICE

Technical Support and E-Warranty Certificate

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Importerad till AUS: SIHAO PTY LTD. 1 ROKEVA STREETEASTWOOD NSW 2122 Australien
Importerad till USA: Sanven Technology Ltd. Suite 250, 9166 Anaheim Place, Rancho Cucamonga, CA 91730
E-CrossStu GmbH