VEVOR[®]

TOUGH TOOLS, HALF PRICE

Technical Support and E-Warranty Certificate www.vevor.com/support

Baby Scale

Model: B18H

We continue to be committed to provide you tools with competitive price. "Save Half", "Half Price" or any other similar expressions used by us only represents an estimate of savings you might benefit from buying certain tools with us compared to the major top brands and does not necessarily mean to cover all categories of tools offered by us. You are kindly reminded to verify carefully when you are placing an order with us if you are actually saving half in comparison with the top major brands.

NEED HELP? CONTACT US!

Have product questions? Need technical support? Please feel free to contact us:

Technical Support and E-Warranty Certificate

www.vevor.com/support

This is the original instruction, please read all manual instructions carefully before operating. VEVOR reserves a clear interpretation of our user manual. The appearance of the product shall be subject to the product you received. Please forgive us that we won't inform you again if there are any technology or software updates on our product.

Thanks for purchasing our baby scale. In order to ensure the accuracy of weighing and its service life, please read this instruction manual carefully before using. This baby scale provides accurate weight and height measurement for baby growth. Suitable for infant and toddlers.

PRODUCT PARAMETER

- High precision sensors
- LED display with white word
- Automatic shutdown
- Unit conversion function
- Tare function
- Hold function
- Overload indication-Err
- Low voltage indication-LO
- HIPS platform
- Large weight measurement range: 30g (0.07lb) - 30kg (66lb)
- Large height measurement range: 0-60cm (23.6 inch)

POWER SAVING MODE

If there is no action within 2 minutes (in a stationary state), the scale will automatically shut down.

BATTERY INSTALLATION

1. Tear off the transparent insulation film on the battery, open the battery cover at the bottom.
2. Please pay attention to the positive and negative pole directions during the installation. Please refer to the “+”、“-” on the battery cover.(See the picture below)
3. Install the battery cover.


Install Scale Guide


Step 1: Install the tray

Step 2: Insert the tray into Scale body

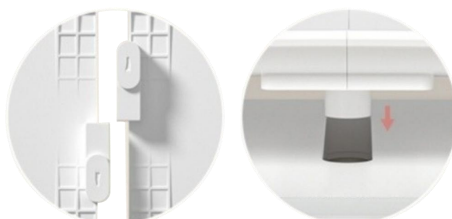Step 3: Twist lock to lock tray

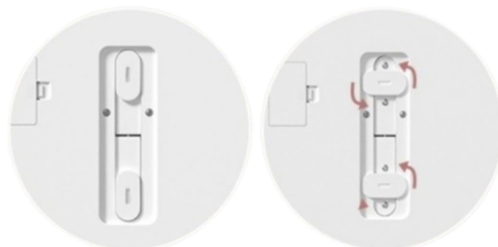
Using Guide

1. Place the scale on a stable, sturdy and level platform.(Do not use it on a shaking/uneven platform.)

2. Short Press “” button to turn it on, after 2s, the display will show “0”;

3. Short press “UNIT” to select the units (kg, oz, lb, lb:oz, g)

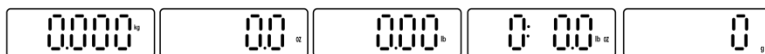4. You can put the Blanket or cotton pad etc on the tray, short press “” to zero out.

5. Put the baby on the scale and press “HOLD” button to lock the weight. (When the display shows “HOLD”, the weighing function is not working unless you press the “HOLD” again to unlock the weight.)

6. After weighing, long press “” to turn off the device or it will automatically turn off after two minutes.

Notice: If you don't use a tray, please let your child stand or sit on the scale steadily.

Notice

1. Choose a suitable floor: It should be on a stable&solid floor or table, and avoid using it in a shaking or uneven place.

2. Be gentle when placing and taking out the baby, and avoid collision or impact on the baby scale.

3. Cleaning: Clean the scale plate frequently to keep it clean . It can be disinfected before use to prevent cross infection.

4. Secure fixation: Make sure that the baby is properly fixed on the scale to prevent it from slipping or rolling.

5. Avoid overload: Do not exceed the maximum load capacity of the baby scale.
6. Battery maintenance: If it is battery-powered, pay attention to replace the battery in time to ensure that the battery is sufficient.
7. Stay away from water sources: Prevent water from splashing on the baby scale to avoid damaging the electronic components.
8. Check the appearance: Check the baby scale for signs of damage or malfunction before use. Do not use it if there is any abnormality.
9. Accurate reading: Record the reading after it stabilizes to avoid reading inaccurate data too early.
10. Storage conditions: When not in use, it should be stored in a suitable environment to avoid extreme environments such as moisture, high temperature or low temperature that may affect its performance.

Safety Information

1. When using, please place the scale on a flat and solid surface, and avoid using it in a shaking or uneven place, which may cause safety problems for infants and children and the scale.
2. Do not use this product in places with high humidity, high altitude, extremely hot or extremely cold.
3. Do not damage the product due to accidental drops, vibrations or other impacts.
4. Please clean it with a soft clean cloth and soapy water in time after use. Do not use chemical solutions such as bleach or alcohol for cleaning.
5. This product is for home use only. This product is not intended for weighing and measuring for commercial purposes.
6. If it is not used for a long time, please remove the battery.
7. If "L0" is displayed, please replace the battery in time.
8. If "Err" overload is displayed, please remove or reduce the object from the scale.

INCLUDES

1. Baby Scale
2. Tray Assembly
3. Instruction Manual
4. 3 AAA Batteries

FCC Information:

CAUTION: Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment!

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- 1) This product may cause harmful interference.
- 2) This product must accept any interference received, including interference that may cause undesired operation.

WARNING: Changes or modifications to this product not expressly approved by the party responsible for compliance could void the user's authority to operate the product.

Note: This product has been tested and found to comply with the limits for a Class B digital device pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This product generates, uses and can radiate radio frequency energy, and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this product does cause harmful interference to radio or television reception, which can be determined by turning the product off and on, the user is encouraged to try to correct the interference by one or more of the following measures.

- Reorient or relocate the receiving antenna.
- Increase the distance between the product and receiver.
- Connect the product to an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for assistance.

CORRECT DISPOSAL:

This product is subject to the provision of European Directive 2012/19/EU. The symbol showing a wheeled bin crossed through indicates that the product requires separate refuse collection in the European Union. This applies to the product and all accessories marked with this symbol. Products marked as such may not be discarded with normal domestic waste, but must be taken to a collection point for recycling electrical and electronic devices.

Manufacturer: Shanghaimuxinmuyeyouxiangongsi
Address: Shuangchenglu 803nong11hao1602A-1609shi, baoshanqu, shanghai 200000 CN.
Imported to AUS: SIHAO PTY LTD. 1 ROKEVA STREETEASTWOOD NSW 2122 Australia
Imported to USA: Sanven Technology Ltd. Suite 250, 9166 Anaheim Place, Rancho Cucamonga, CA 91730

UK	REP
----	-----

YH CONSULTING LIMITED.
C/O YH Consulting Limited Office 147,
Centurion House, London Road,
Staines-upon-Thames, Surrey, TW18 4AX

EC	REP
----	-----

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Bluetooth Baby Scale Instruction

Free BabyFit APP

Search "BabyFit" in the application market.

get on Google Play

get on Google Play

You can scan the QR code to download the app, or search APP named "BabyFit" in the app store to download it.


Please read the instruction before using this item
BabyFit downloading guide

You can scan the QR code to download the app, or search APP named “BabyFit “in the app store to download it.


After download the app , open the app and follow the prompts to register information.

Steps guide

Step 1: Create a new account (you can use your email /Facebook /Wechat to creative your babyFit account)


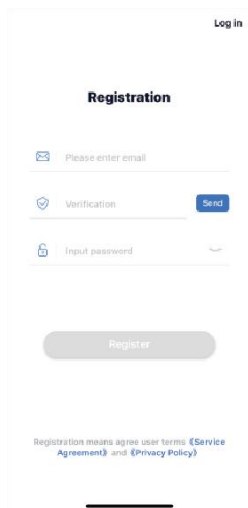
The screenshot shows the 'Registration' screen of the BabyFit app. At the top right is a 'Log in' link. The main heading is 'Registration'. Below it are three input fields: 'Please enter email' with an envelope icon, 'Verification' with a shield icon and a blue 'Send' button, and 'Input password' with a lock icon and a show/hide toggle. A large grey 'Register' button is at the bottom. At the very bottom, there is a line of text: 'Registration means agree user terms 《Service Agreement》 and 《Privacy Policy》'.

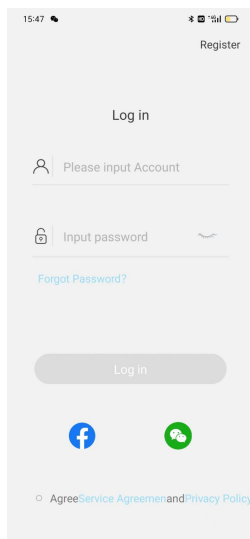The screenshot shows the 'Log in' screen of the BabyFit app. At the top right is a 'Register' link. The main heading is 'Log in'. Below it are two input fields: 'Please input Account' with a person icon and 'Input password' with a lock icon and a show/hide toggle. A blue link 'Forgot Password?' is below the password field. A large grey 'Log in' button is at the bottom. Below the button are two social media icons: Facebook and WeChat. At the very bottom, there is a line of text: '○ Agree Service Agreement and Privacy Policy'.

Step 2: after registration, turn on your phone's Bluetooth and turn on the baby scale , then click “add device”. When device icon appears on the screen and click the icon to pair successfully.

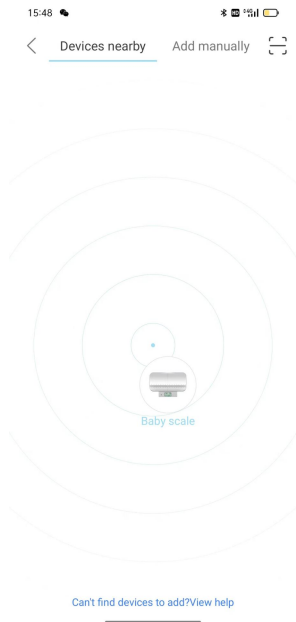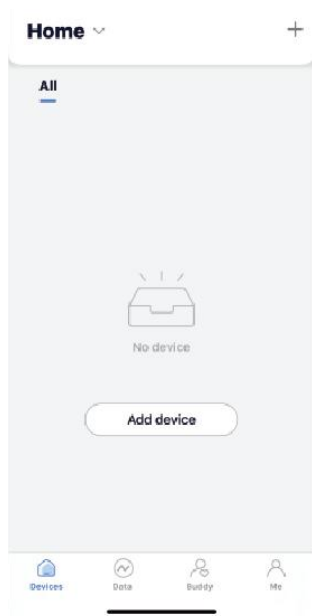Step 3: Set the device name and using location

15:48

Bluetooth, NO, 4G, and battery status icons

Device name

Baby scale >

Choose device place

Living room

Bedroom

Kitchen

Dining room

Study room

Bathroom

Start to use

Step 4: setting weight and height unit, choose the unit you need

<

Unit setting

Weight unit

kg ✓

lb:oz

oz

lb


Height unit

cm ✓

inch

OK


Step 5: Enter a nickname for your baby and select the relationship (son, daughter or other) and select the date of birth and gender (male/female)

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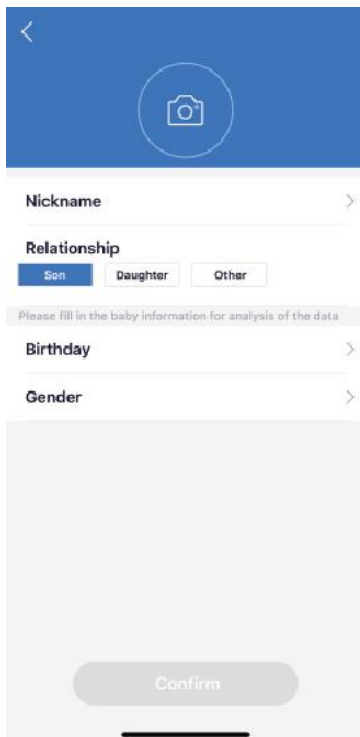- Nickname:** A text input field with a chevron icon on the right.
- Relationship:** Three buttons labeled "Son", "Daughter", and "Other". The "Son" button is highlighted in blue.
- Instructional text:** A light gray bar containing the text "Please fill in the baby information for analysis of the data".
- Birthday:** A text input field with a chevron icon on the right.
- Gender:** A text input field with a chevron icon on the right.
- Confirm button:** A large, light gray rounded rectangle at the bottom with the word "Confirm" in the center.


Step 6: Once the settings are complete, you can weigh your baby. Weigh your baby and click "Record" to save the data. You can manually enter weight, height and head circumference (HG)

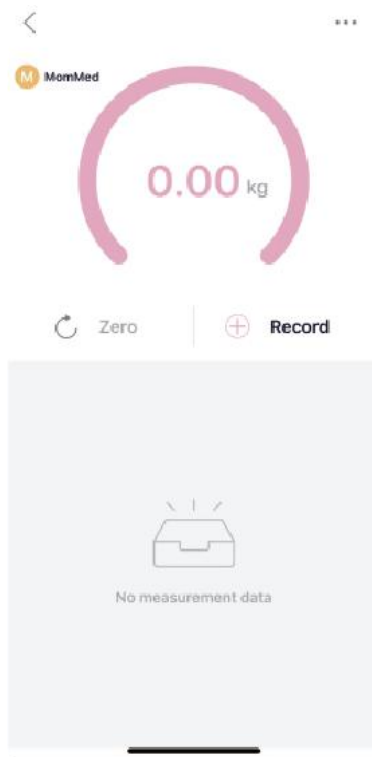Step 7: If you want to view records, click to enter your previous records to view summary information. You can view records of weight, height, head circumference (HC) and body mass index (BMI).

15:49					
<	U29980				
Record	Weight	Height	HC	BMI	
2024-06-07 2Y-11M					
Weight	2.260 kg	Height	20.0 cm		
HC	20.0 cm	BMI	56.5		
Your baby's physical development has not reached the normal range!					
2024-06-03 2Y-11M					
Weight	0.200 kg	Height	55.70 cm		
HC	37.50 cm	BMI	0.6		
Your baby's physical development has not reached the normal range!					
2024-05-23 2Y-11M					
Weight	0.51 kg	Height	55.70 cm		
HC	37.50 cm	BMI	1.6		
Your baby's physical development has not reached the normal range!					
2024-05-23 2Y-11M					

Step 8: For the growth curve, you can view the data points on the chart.

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
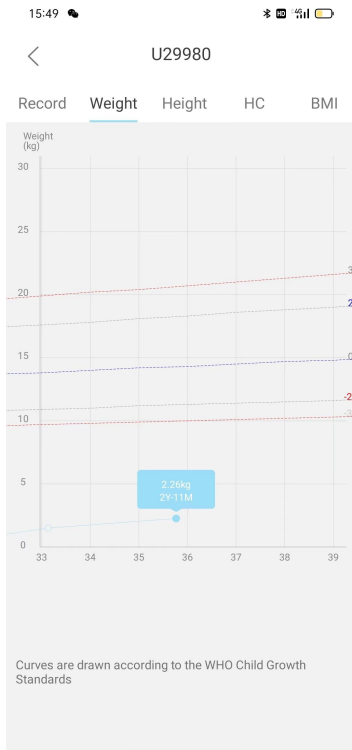
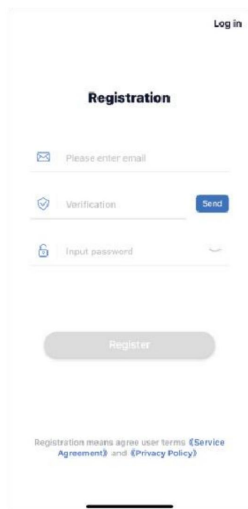
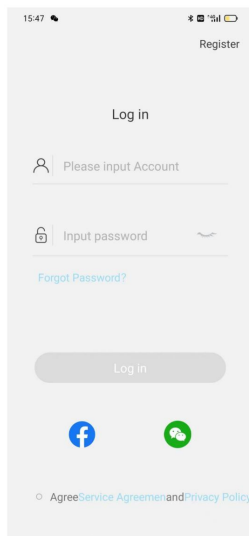
The screenshot shows the 'Login' screen of the BabyFit app. At the top right is a 'Register' link. The main heading is 'Log in'. Below it are two input fields: 'Please input Account' with a person icon and 'Input password' with a lock icon and a show/hide toggle. A blue link 'Forgot Password?' is below the password field. A large grey 'Log in' button is at the bottom. Below the button are two social media icons: Facebook and WeChat. At the very bottom, there is a line of text: 'Agree Service Agreement and Privacy Policy'.

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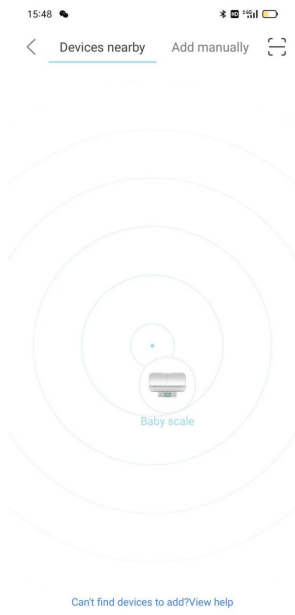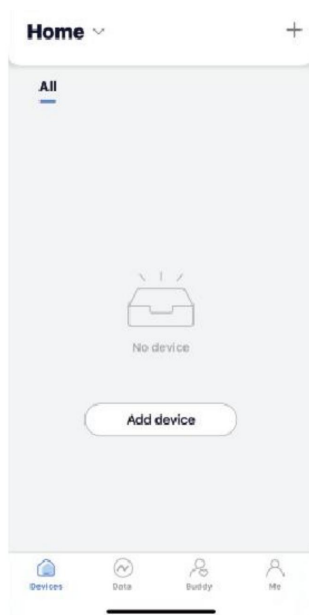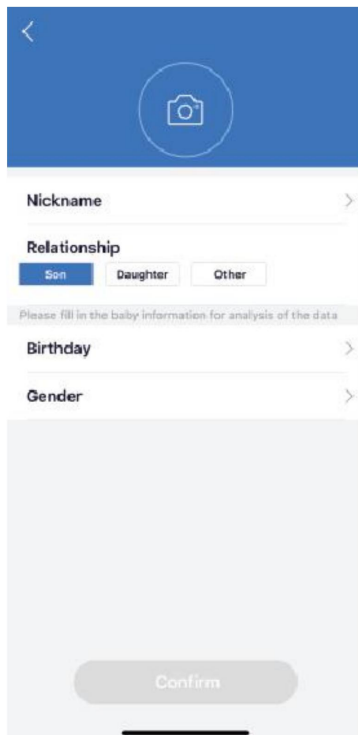
The screenshot shows a mobile application interface for entering baby information. At the top is a blue header with a back arrow on the left and a camera icon in a circle in the center. Below the header, the form has several sections: a 'Nickname' text input field with a right arrow; a 'Relationship' section with three buttons: 'Son' (highlighted in blue), 'Daughter', and 'Other'; a light gray instruction bar that says 'Please fill in the baby information for analysis of the data'; a 'Birthday' text input field with a right arrow; and a 'Gender' text input field with a right arrow. At the bottom of the form is a large, light gray 'Confirm' button. The entire form is set against a light gray background.

<

Camera icon

Nickname >

Relationship

Son Daughter Other


Please fill in the baby information for analysis of the data

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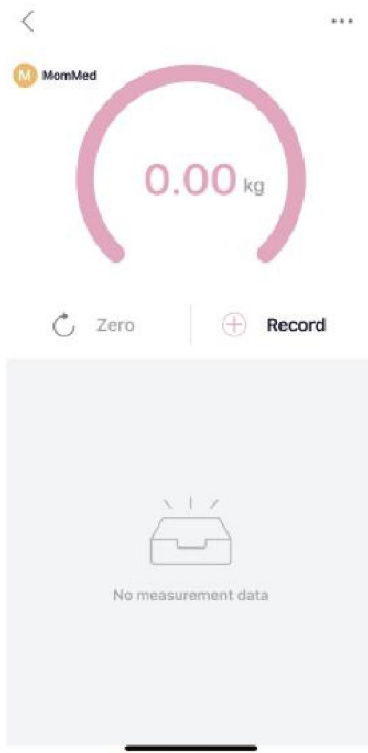15:49		* 4G 5G	
<		U29980	
Record	Weight	Height	HC BMI
2024-06-07 2Y-11M			
Weight	2.260 kg	Height	20.0 cm
HC	20.0 cm	BMI	56.5
Your baby's physical development has not reached the normal range!			
2024-06-03 2Y-11M			
Weight	0.200 kg	Height	55.70 cm
HC	37.50 cm	BMI	0.6
Your baby's physical development has not reached the normal range!			
2024-05-23 2Y-11M			
Weight	0.51 kg	Height	55.70 cm
HC	37.50 cm	BMI	1.6
Your baby's physical development has not reached the normal range!			
2024-05-23 2Y-11M			

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
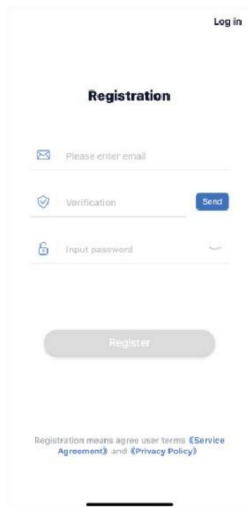
The screenshot shows the 'Registration' screen of the BabyFit app. At the top right is a 'Log in' link. The main heading is 'Registration'. Below it are three input fields: 'Please enter email' with an envelope icon, 'Verification' with a shield icon and a blue 'Send' button, and 'Input password' with a lock icon and a show/hide toggle. A large grey 'Register' button is at the bottom. At the very bottom, there is a line of text: 'Registration means agree user terms (Service Agreement) and (Privacy Policy)'.


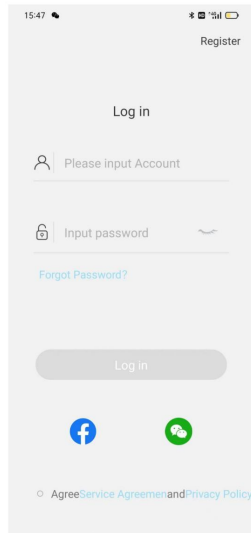
The screenshot shows the 'Log in' screen of the BabyFit app. At the top right is a 'Register' link. The main heading is 'Log in'. Below it are two input fields: 'Please input Account' with a person icon and 'Input password' with a lock icon and a show/hide toggle. A blue link 'Forgot Password?' is below the password field. A large grey 'Log in' button is at the bottom. Below the button are two circular icons for Facebook and WeChat. At the very bottom, there is a line of text: 'Agree Service Agreement and Privacy Policy'.

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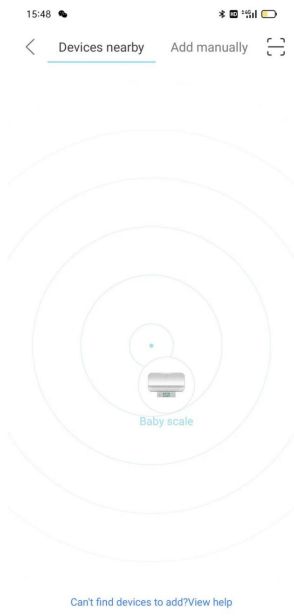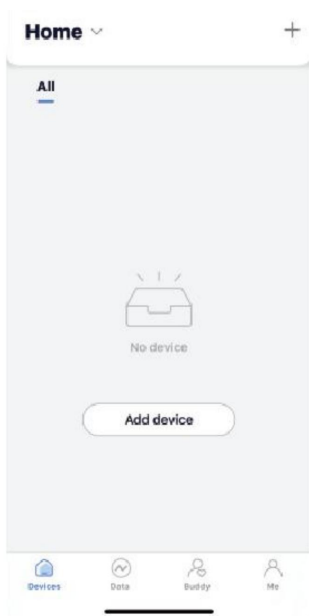Device name Baby scale >

Choose device place

Living room Bedroom Kitchen Dining room

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
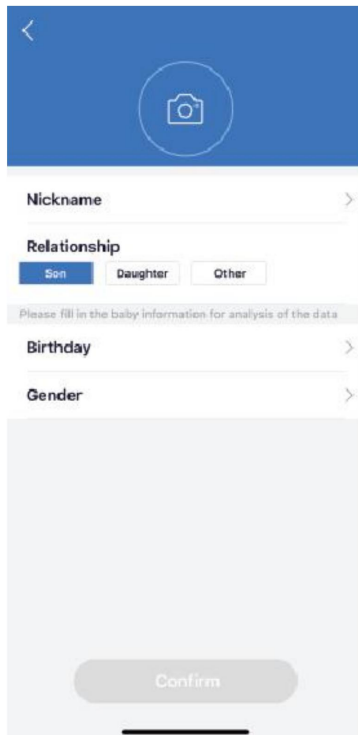
The screenshot shows a mobile application interface for entering baby information. At the top is a blue header with a back arrow on the left and a camera icon in a circle in the center. Below the header, the form contains the following fields and options:

- Nickname:** A text input field with a right-pointing chevron.
- Relationship:** A section with three buttons: "Son" (highlighted in blue), "Daughter", and "Other".
- Instruction:** A light gray bar with the text "Please fill in the baby information for analysis of the data".
- Birthday:** A text input field with a right-pointing chevron.
- Gender:** A text input field with a right-pointing chevron.
- Confirm:** A large, light gray rounded button at the bottom of the form.

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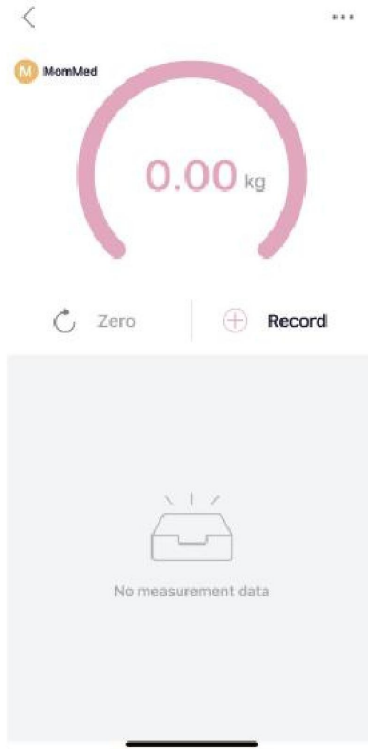
15:49		* 5G 4G	
<		U29980	
Record	Weight	Height	HC BMI
2024-06-07 2Y-11M			
Weight	2.260 kg	Height	20.0 cm
HC	20.0 cm	BMI	56.5
Your baby's physical development has not reached the normal range!			
2024-06-03 2Y-11M			
Weight	0.200 kg	Height	55.70 cm
HC	37.50 cm	BMI	0.6
Your baby's physical development has not reached the normal range!			
2024-05-23 2Y-11M			
Weight	0.51 kg	Height	55.70 cm
HC	37.50 cm	BMI	1.6
Your baby's physical development has not reached the normal range!			
2024-05-23 2Y-11M			

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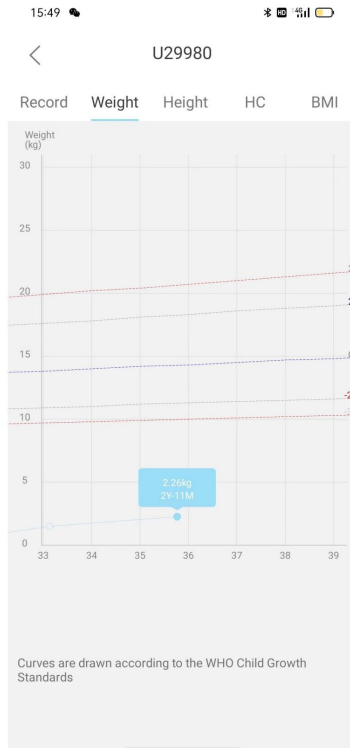
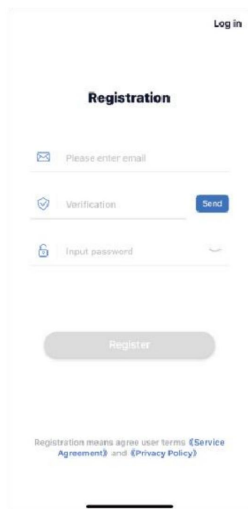
The screenshot shows the 'Registration' screen of the BabyFit app. At the top right is a 'Log in' link. The main title is 'Registration'. Below it are three input fields: 'Please enter email' with an envelope icon, 'Verification' with a shield icon and a blue 'Send' button, and 'Input password' with a lock icon and a show/hide toggle. A large grey 'Register' button is at the bottom. At the very bottom, there is a line of text: 'Registration means agree user terms (Service Agreement) and (Privacy Policy)'.


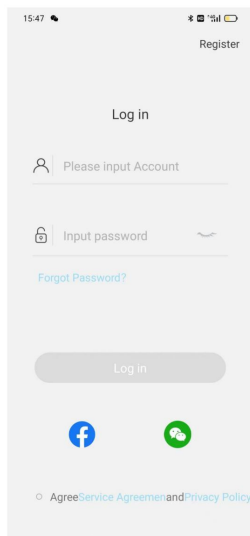
The screenshot shows the 'Log in' screen of the BabyFit app. At the top right is a 'Register' link. The main title is 'Log in'. Below it are two input fields: 'Please input Account' with a person icon and 'Input password' with a lock icon and a show/hide toggle. A blue link 'Forgot Password?' is below the password field. A large grey 'Log in' button is at the bottom. Below the button are two social media icons: Facebook and WeChat. At the very bottom, there is a line of text: 'Agree Service Agreement and Privacy Policy'.

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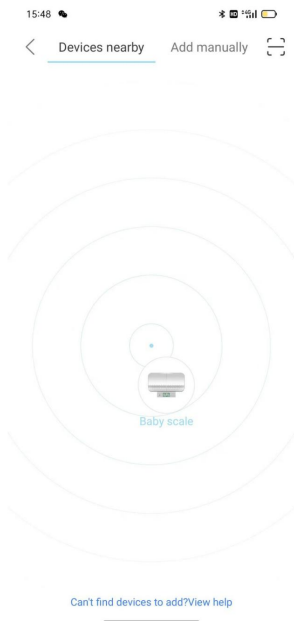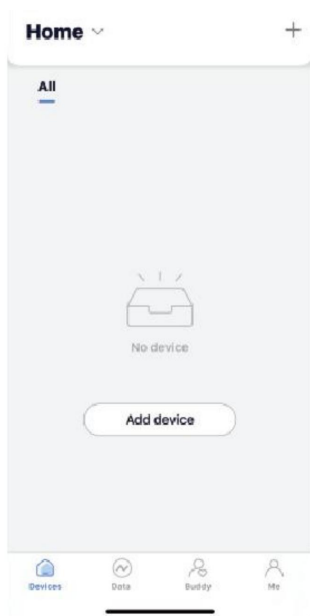The image shows a mobile application screen titled "Unit setting". At the top left is a back arrow icon. The screen is divided into two main sections: "Weight unit" and "Height unit".

Weight unit section:


- The "kg" option is selected, indicated by a checkmark icon to its right.
- Below "kg" are three unselected options: "lb:oz", "oz", and "lb", each with a horizontal line for selection.

Height unit section:

- The "cm" option is selected, indicated by a checkmark icon to its right.
- Below "cm" is one unselected option: "inch", with a horizontal line for selection.

At the bottom of the screen is a large, rounded rectangular button labeled "OK". Below the "OK" button is a thick black horizontal line, likely representing the home indicator bar on an iPhone.

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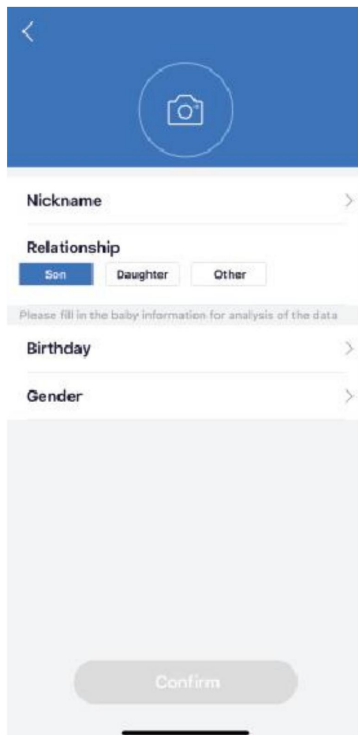
The screenshot shows a mobile application interface for entering baby information. At the top is a blue header with a back arrow on the left and a camera icon in a circle in the center. Below the header is a form with the following sections:

- Nickname:** A text input field with a right-pointing chevron.
- Relationship:** Three buttons: "Son" (highlighted in blue), "Daughter", and "Other".
- A light gray separator bar with the text: "Please fill in the baby information for analysis of the data".
- Birthday:** A text input field with a right-pointing chevron.
- Gender:** A text input field with a right-pointing chevron.

At the bottom of the form is a large, light gray rectangular area, and at the very bottom is a rounded "Confirm" button. A black horizontal line is visible at the bottom edge of the screen, indicating a mobile device interface.

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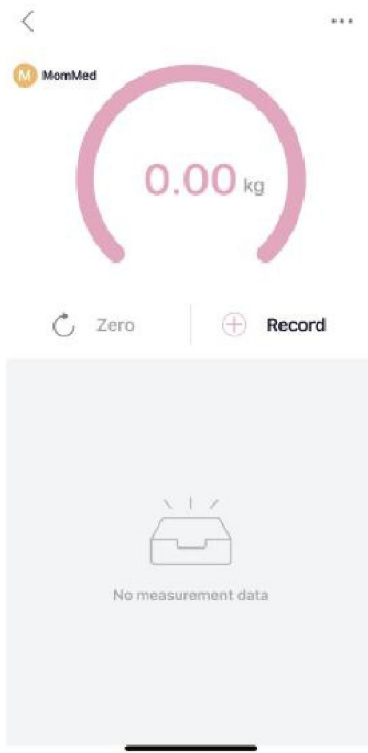BMI

56.5

Your baby's physical development has not reached the normal range!

2024-06-03 2Y-11M

Weight

0.200

kg

Height

55.70

cm

HC

37.50

cm

BMI

0.6

Your baby's physical development has not reached the normal range!

2024-05-23 2Y-11M

Weight

0.51

kg

Height

55.70

cm

HC

37.50

cm


BMI

1.6

Your baby's physical development has not reached the normal range!

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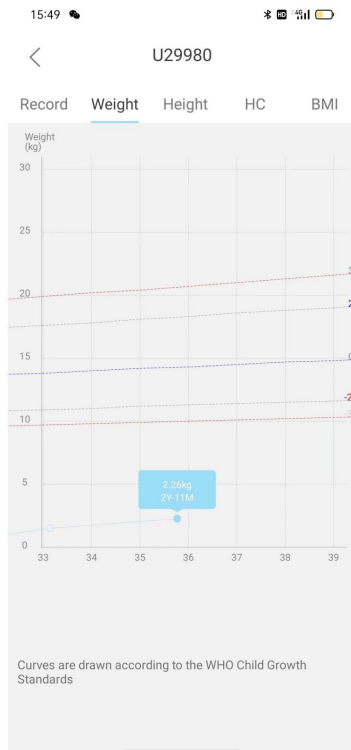
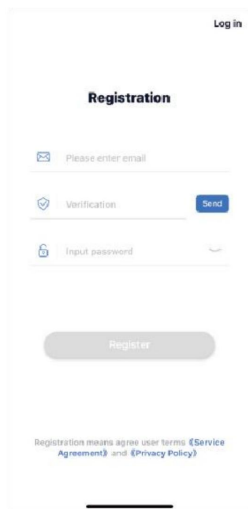
This screenshot shows the registration interface of the BabyFit app. At the top right is a "Log in" link. The main heading is "Registration". Below it are three input fields: "Please enter email" with an envelope icon, "Verification" with a shield icon and a blue "Send" button, and "Input password" with a lock icon and a show/hide toggle. A large grey "Register" button is at the bottom. At the very bottom, there is a line of text: "Registration means agree user terms (Service Agreement) and (Privacy Policy)".


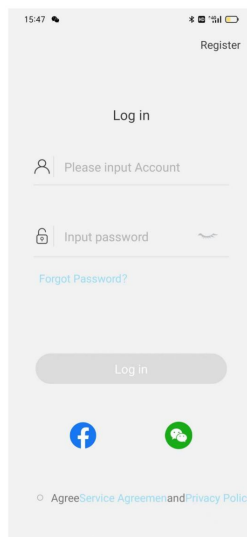
This screenshot shows the login interface of the BabyFit app. At the top right is a "Register" link. The main heading is "Log in". Below it are two input fields: "Please input Account" with a person icon and "Input password" with a lock icon and a show/hide toggle. A blue link "Forgot Password?" is below the password field. A large grey "Log in" button is at the bottom. Below the button are two social media icons: Facebook and WeChat. At the very bottom, there is a line of text: "Agree Service Agreement and Privacy Policy".

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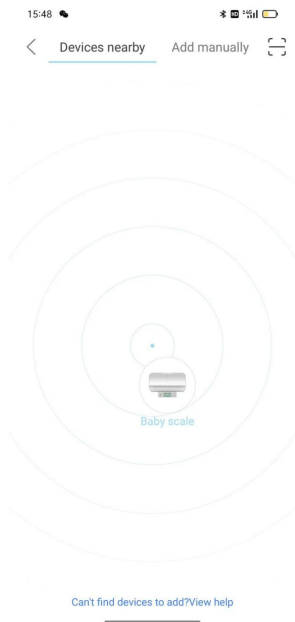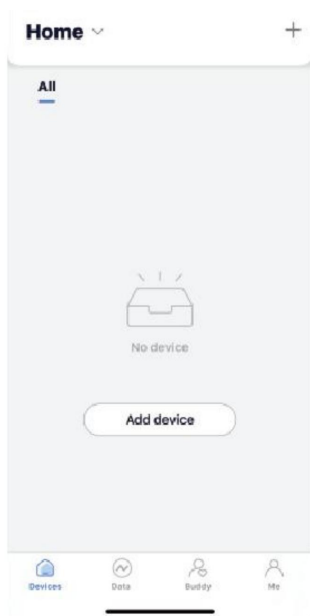
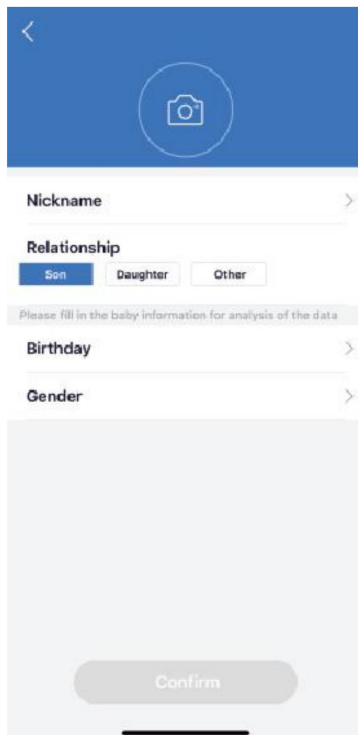
The screenshot shows a mobile application interface for adding a baby's information. At the top, there is a blue header with a back arrow on the left and a camera icon in the center, indicating a profile picture. Below the header, the form contains the following fields and options:

- Nickname:** A text input field with a chevron icon on the right.
- Relationship:** A section with three buttons: "Son" (selected, highlighted in blue), "Daughter", and "Other".
- Instruction:** A light gray bar with the text "Please fill in the baby information for analysis of the data".
- Birthday:** A text input field with a chevron icon on the right.
- Gender:** A text input field with a chevron icon on the right.
- Confirm:** A large, rounded gray button at the bottom of the form.

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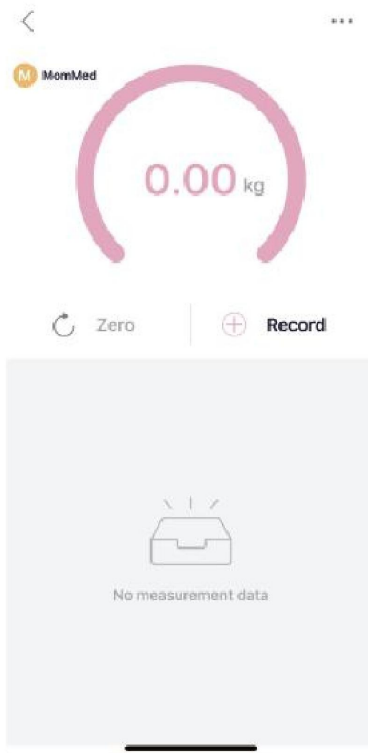
15:49				
<	U29980			
Record	Weight	Height	HC	BMI
2024-06-07 2Y-11M				
Weight	2.260 kg	Height	20.0 cm	
HC	20.0 cm	BMI	56.5	
Your baby's physical development has not reached the normal range!				
2024-06-03 2Y-11M				
Weight	0.200 kg	Height	55.70 cm	
HC	37.50 cm	BMI	0.6	
Your baby's physical development has not reached the normal range!				
2024-05-23 2Y-11M				
Weight	0.51 kg	Height	55.70 cm	
HC	37.50 cm	BMI	1.6	
Your baby's physical development has not reached the normal range!				
2024-05-23 2Y-11M				

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

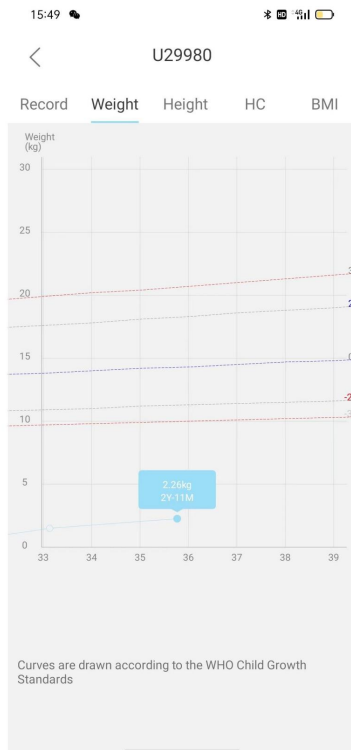
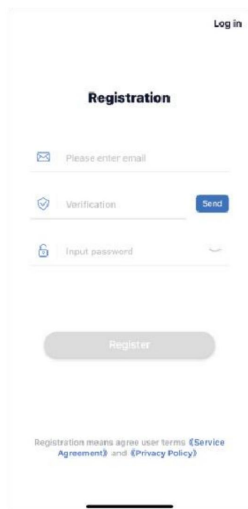
The screenshot shows the 'Registration' screen of the BabyFit app. At the top right is a 'Log in' link. The main heading is 'Registration'. Below it are three input fields: 'Please enter email' with an envelope icon, 'Verification' with a shield icon and a blue 'Send' button, and 'Input password' with a lock icon and a show/hide toggle. A large grey 'Register' button is at the bottom. At the very bottom, there is a line of text: 'Registration means agree user terms (Service Agreement) and (Privacy Policy)'.

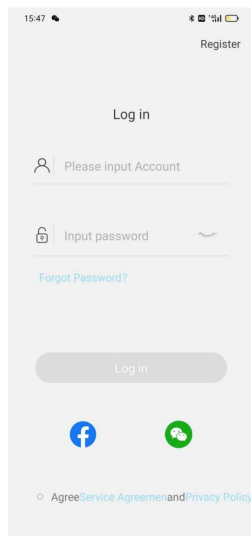The screenshot shows the 'Log in' screen of the BabyFit app. At the top right is a 'Register' link. The main heading is 'Log in'. Below it are two input fields: 'Please input Account' with a person icon and 'Input password' with a lock icon and a show/hide toggle. A blue link 'Forgot Password?' is below the password field. A large grey 'Log in' button is at the bottom. Below the button are two social media icons: Facebook and WeChat. At the very bottom, there is a line of text: 'Agree Service Agreement and Privacy Policy'.

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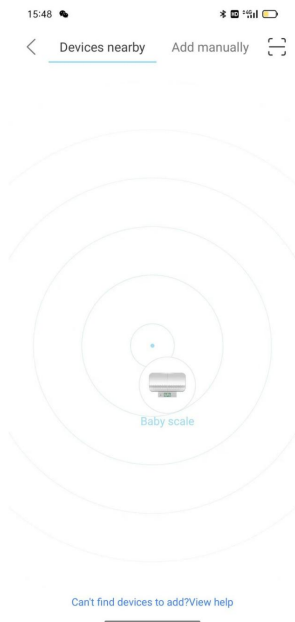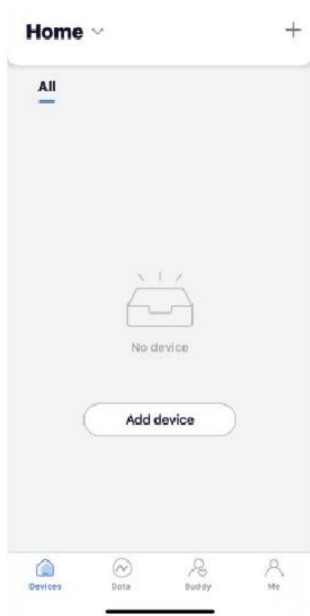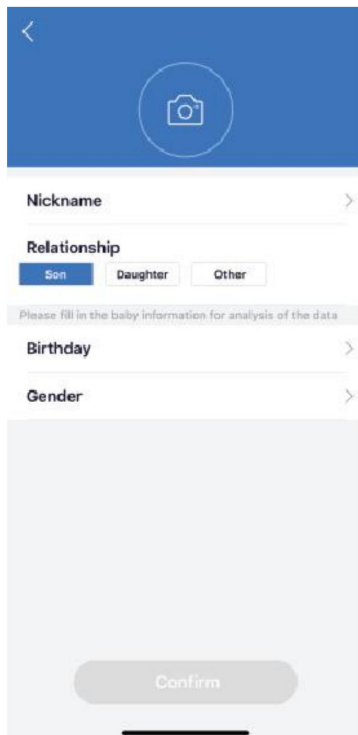
The screenshot shows a mobile application interface for adding a baby. At the top is a blue header with a back arrow on the left and a camera icon in a circle in the center. Below the header, the form contains the following fields and options:

- Nickname**: A text input field with a right-pointing chevron.
- Relationship**: Three buttons labeled "Son", "Daughter", and "Other". The "Son" button is currently selected.
- A grey separator bar with the text: "Please fill in the baby information for analysis of the data".
- Birthday**: A text input field with a right-pointing chevron.
- Gender**: A text input field with a right-pointing chevron.
- A large, light grey rectangular area below the input fields.
- A rounded "Confirm" button at the bottom of the form.

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

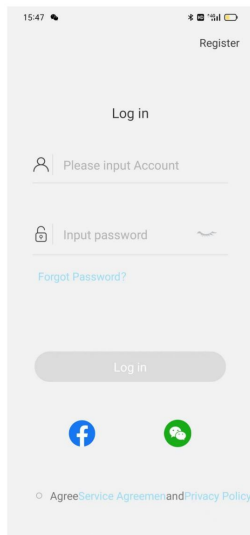
This screenshot shows the login interface of the BabyFit app. At the top right is a "Register" link. The main heading is "Log in". Below it are two input fields: "Please input Account" with a person icon and "Input password" with a lock icon and a show/hide toggle. A blue link "Forgot Password?" is below the password field. A large grey "Log in" button is at the bottom. Below the button are two social media icons: Facebook and WeChat. At the very bottom, there is a line of text: "○ Agree Service Agreement and Privacy Policy".

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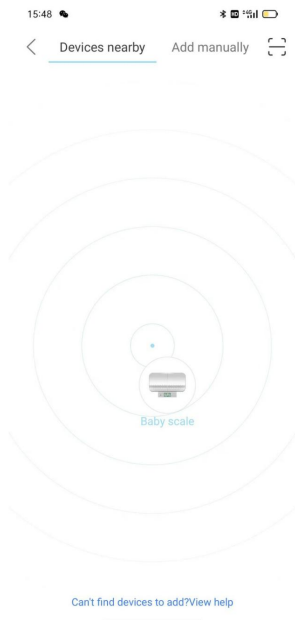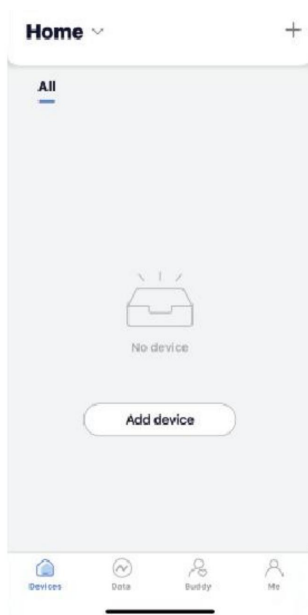
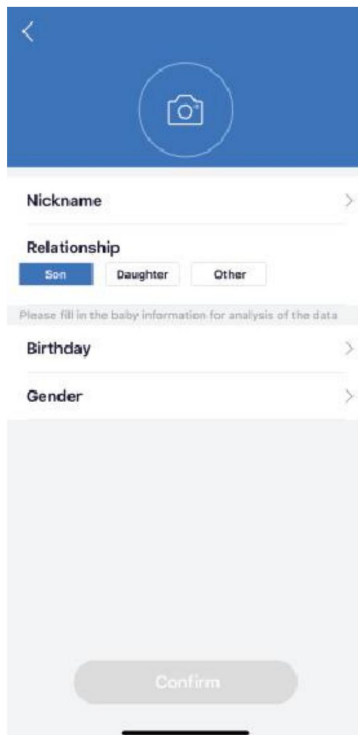
The screenshot shows a mobile application interface for entering baby information. At the top is a blue header with a back arrow on the left and a camera icon in a circle in the center. Below the header, there is a form with the following sections:

- Nickname:** A text input field with a right-pointing chevron.
- Relationship:** Three buttons: "Son" (highlighted in blue), "Daughter", and "Other".
- A grey separator bar with the text: "Please fill in the baby information for analysis of the data".
- Birthday:** A text input field with a right-pointing chevron.
- Gender:** A text input field with a right-pointing chevron.

At the bottom of the form is a large, light grey rectangular area, and at the very bottom is a rounded "Confirm" button.

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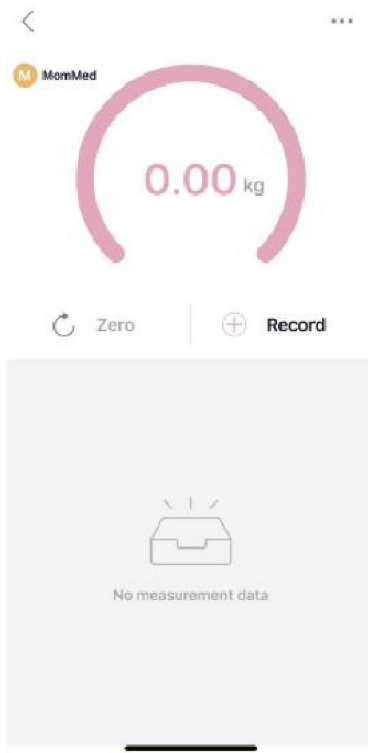15:49		* 4G LTE	
<		U29980	
Record	Weight	Height	HC BMI
2024-06-07 2Y-11M			
Weight	2.260 kg	Height	20.0 cm
HC	20.0 cm	BMI	56.5
Your baby's physical development has not reached the normal range!			
2024-06-03 2Y-11M			
Weight	0.200 kg	Height	55.70 cm
HC	37.50 cm	BMI	0.6
Your baby's physical development has not reached the normal range!			
2024-05-23 2Y-11M			
Weight	0.51 kg	Height	55.70 cm
HC	37.50 cm	BMI	1.6
Your baby's physical development has not reached the normal range!			
2024-05-23 2Y-11M			

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
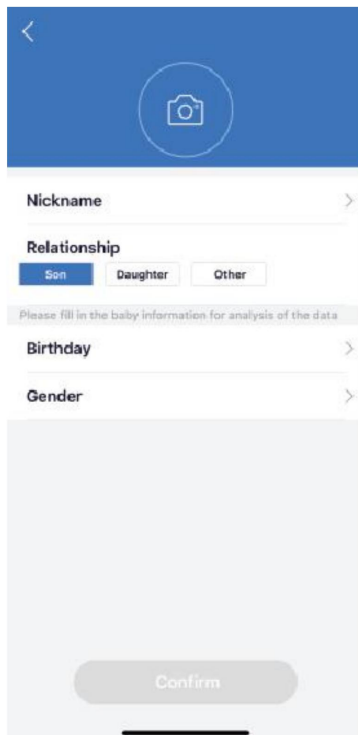
The screenshot shows a mobile application interface for adding a new baby. At the top is a blue header with a back arrow on the left and a camera icon in a circle in the center. Below the header, the form contains the following fields and options:

- Nickname:** A text input field with a chevron icon on the right.
- Relationship:** A section with three buttons: "Son" (highlighted in blue), "Daughter", and "Other".
- Instruction:** A light gray bar with the text "Please fill in the baby information for analysis of the data".
- Birthday:** A text input field with a chevron icon on the right.
- Gender:** A text input field with a chevron icon on the right.
- Confirm:** A large, rounded gray button at the bottom of the form.

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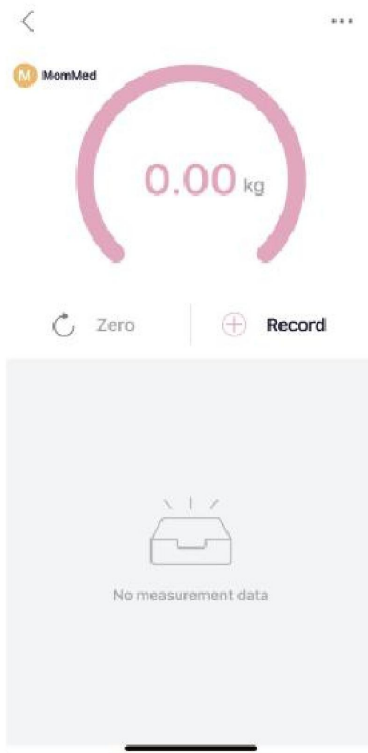15:49		U29980	
<		✎	
Record	Weight	Height	HC BMI
2024-06-07 2Y-11M			
Weight	2.260 kg	Height	20.0 cm
HC	20.0 cm	BMI	56.5
Your baby's physical development has not reached the normal range!			
2024-06-03 2Y-11M			
Weight	0.200 kg	Height	55.70 cm
HC	37.50 cm	BMI	0.6
Your baby's physical development has not reached the normal range!			
2024-05-23 2Y-11M			
Weight	0.51 kg	Height	55.70 cm
HC	37.50 cm	BMI	1.6
Your baby's physical development has not reached the normal range!			
2024-05-23 2Y-11M			

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