VEVOR[®]

TOUGH TOOLS, HALF PRICE

Technical Support and E-Warranty Certificate www.vevor.com/support

Baby Scale

Model: B18H

We continue to be committed to provide you tools with competitive price. "Save Half", "Half Price" or any other similar expressions used by us only represents an estimate of savings you might benefit from buying certain tools with us compared to the major top brands and does not necessarily mean to cover all categories of tools offered by us. You are kindly reminded to verify carefully when you are placing an order with us if you are actually saving half in comparison with the top major brands.

VEVOR[®]

TOUGH TOOLS, HALF PRICE

Baby Scale Instruction Manual

NEED HELP? CONTACT US!

Have product questions? Need technical support? Please feel free to contact us:

Technical Support and E-Warranty Certificate

www.vevor.com/support

This is the original instruction, please read all manual instructions carefully before operating. VEVOR reserves a clear interpretation of our user manual. The appearance of the product shall be subject to the product you received. Please forgive us that we won't inform you again if there are any technology or software updates on our product.

Thanks for purchasing our baby scale. In order to ensure the accuracy of weighing and its service life, please read this instruction manual carefully before using. This baby scale provides accurate weight and height measurement for baby growth. Suitable for infant and toddlers.

PRODUCT PARAMETER


- High precision sensors
- LED display with white word
- Automatic shutdown
- Unit conversion function
- Tare function
- Hold function
- Overload indication-Err
- Low voltage indication-LO
- HIPS platform
- Large weight measurement range: 30g (0.07lb) - 30kg (66lb)
- Large height measurement range: 0-60cm (23.6 inch)

POWER SAVING MODE

If there is no action within 2 minutes (in a stationary state), the scale will automatically shut down.

BATTERY INSTALLATION

1. Tear off the transparent insulation film on the battery, open the battery cover at the bottom.
2. Please pay attention to the positive and negative pole directions during the installation. Please refer to the “+”、“-” on the battery cover.(See the picture below)
3. Install the battery cover.

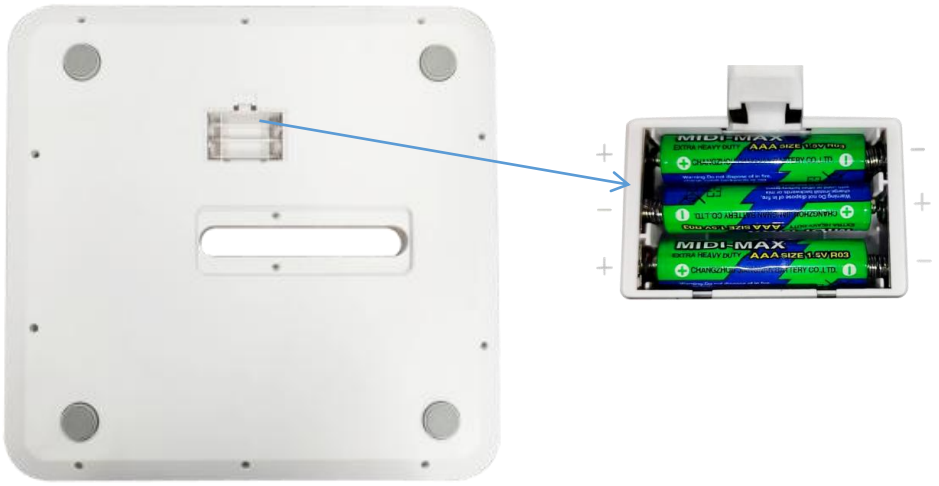
Install Scale Guide


Step 1: Install the tray

Step 2: Insert the tray into Scale body

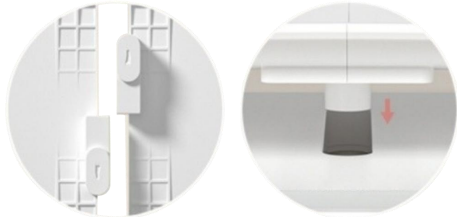Step 3: Twist lock to lock tray

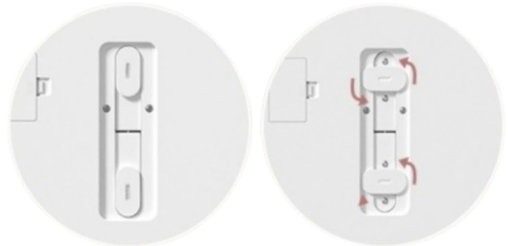
Using Guide

1. Place the scale on a stable, sturdy and level platform.(Do not use it on a shaking/uneven platform.)

2. Short Press “” button to turn it on, after 2s, the display will show “0”;

3. Short press “UNIT” to select the units (kg, oz, lb, lb:oz, g)

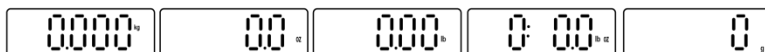4. You can put the Blanket or cotton pad etc on the tray, short press “” to zero out.

5. Put the baby on the scale and press “HOLD” button to lock the weight. (When the display shows “HOLD”, the weighing function is not working unless you press the “HOLD” again to unlock the weight.)

6. After weighing, long press “” to turn off the device or it will automatically turn off after two minutes.

Notice: If you don't use a tray, please let your child stand or sit on the scale steadily.

Notice

1. Choose a suitable floor: It should be on a stable&solid floor or table, and avoid using it in a shaking or uneven place.

2. Be gentle when placing and taking out the baby, and avoid collision or impact on the baby scale.

3. Cleaning: Clean the scale plate frequently to keep it clean . It can be disinfected before use to prevent cross infection.

4. Secure fixation: Make sure that the baby is properly fixed on the scale to prevent it from slipping or rolling.

5. Avoid overload: Do not exceed the maximum load capacity of the baby scale.
6. Battery maintenance: If it is battery-powered, pay attention to replace the battery in time to ensure that the battery is sufficient.
7. Stay away from water sources: Prevent water from splashing on the baby scale to avoid damaging the electronic components.
8. Check the appearance: Check the baby scale for signs of damage or malfunction before use. Do not use it if there is any abnormality.
9. Accurate reading: Record the reading after it stabilizes to avoid reading inaccurate data too early.
10. Storage conditions: When not in use, it should be stored in a suitable environment to avoid extreme environments such as moisture, high temperature or low temperature that may affect its performance.

Safety Information

1. When using, please place the scale on a flat and solid surface, and avoid using it in a shaking or uneven place, which may cause safety problems for infants and children and the scale.
2. Do not use this product in places with high humidity, high altitude, extremely hot or extremely cold.
3. Do not damage the product due to accidental drops, vibrations or other impacts.
4. Please clean it with a soft clean cloth and soapy water in time after use. Do not use chemical solutions such as bleach or alcohol for cleaning.
5. This product is for home use only. This product is not intended for weighing and measuring for commercial purposes.
6. If it is not used for a long time, please remove the battery.
7. If "L0" is displayed, please replace the battery in time.
8. If "Err" overload is displayed, please remove or reduce the object from the scale.

INCLUDES

1. Baby Scale
2. Tray Assembly
3. Instruction Manual
4. 3 AAA Batteries

FCC Information:

CAUTION: Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment!

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- 1) This product may cause harmful interference.
- 2) This product must accept any interference received, including interference that may cause undesired operation.

WARNING: Changes or modifications to this product not expressly approved by the party responsible for compliance could void the user's authority to operate the product.

Note: This product has been tested and found to comply with the limits for a Class B digital device pursuant to Part 15 of the FCC Rules, These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This product generates, uses and can radiate radio frequency energy, and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this product does cause harmful interference to radio or television reception, which can be determined by turning the product off and on, the user is encouraged to try to correct the interference by one or more of the following measures.

- Reorient or relocate the receiving antenna.
- Increase the distance between the product and receiver.
- Connect the product to an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for assistance.

CORRECT DISPOSAL:

This product is subject to the provision of European Directive 2012/19/EU. The symbol showing a wheeled bin crossed through indicates that the product requires separate refuse collection in the European Union. This applies to the product and all accessories marked with this symbol. Products marked as such may not be discarded with normal domestic waste, but must be taken to a collection point for recycling electrical and electronic devices.

Manufacturer: Shanghaimuxinmuyeyouxiangongsi

Address: Shuangchenglu 803nong11hao1602A-1609shi, baoshanqu, shanghai 200000 CN.

Imported to AUS: SIHAO PTY LTD. 1 ROKEVA STREETEASTWOOD NSW 2122 Australia

Imported to USA: Sanven Technology Ltd. Suite 250, 9166 Anaheim Place, Rancho Cucamonga, CA 91730

YH CONSULTING LIMITED.
C/O YH Consulting Limited Office 147,
Centurion House, London Road,
Staines-upon-Thames, Surrey, TW18 4AX

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Importerad till AUS: SIHAO PTY LTD. 1 ROKEVA STREETEASTWOOD NSW 2122 Australien
Importerad till USA: Sanven Technology Ltd. Suite 250, 9166 Anaheim Place, Rancho Cucamonga, CA 91730
E-CrossStu GmbH