740mm

210mm

**VEVOR**  
TOUGH TOOLS, HALF PRICE

**Fall Mats**  
Instruction Manual

**NEED HELP? CONTACT US!**  
Have product questions? Need technical support? Please feel free to contact us:  
**Technical Support and E-Warranty Certificate**  
[www.vevor.com/support](http://www.vevor.com/support)

This is the original instruction, please read all manual instructions carefully before operating. VEVOR reserves a clear interpretation of our user manual. The appearance of the product shall be subject to the product you received. Please forgive us that we won't inform you again if there are any technology or software updates on our product.

**Warning:** To reduce the risk of injury, the user must read the instructions manual carefully.

**FEATURES**

- Thank you for using Our company's Wireless off Bed zoning pre-alarm system. Please read the user manual carefully for your first use and keep it for future reference.
- The system unit is small in size and easy to use, with multiple functions such as sound and light warning/vibration reminder /SOS emergency help, bed alarm/near bed warning, etc., which can be applied to the elderly or people with mobility difficulties in homes, nursing homes, hospitals and other places. Wheelchair assistance.
- System unit AA\*2 battery/pressure pad sensor CR2477\*2 battery/Emergency button CR2032\*1 Battery operation, low power consumption, standby time up to one year, with low power reminder function.
- System alarm main unit : High volume + vibration / low volume + vibration / only vibration ,multiple modes are optional, press the MODE key to set the mode you need.

**COMPONENTS**

**INSTALL AND USE**

**For 636**  
90cm x 25cm (for bed)

**Warning one side zone**  
**Warning area on both sides**  
**Pressure pad center**  
Get out of bed alarm sensing area

- Start the system main unit. Transmit sensor, Call button battery cover, load the battery according to the battery +/- polarity mark. (Note: AA\*2 Pcs alkaline battery is suitable for the alarm main unit, AAA \* 2PCS for 636S small transmitter, CR2477\*2 battery is suitable for the pressure pad sensor, CR2032\*1 battery is suitable for the Call button), the main unit power ON/OFF switch is on right side, Set transmit sensor Single / Double and N.W switch option as required.
- System main unit and sensors. The Call button (Optional) has been paired in the factory, the sensor pressure pad is placed flat on the bed or wheelchair, the small bag equipped with the transmitter is hung on the edge of the bed or wheelchair, and the optimal position of the pressure pad is adjusted according to the position of the human body lying or sitting (suggest that the pressure pad be placed flat on the shoulders or under the hip of the human body). Set the N.W/Single/Double switch option of the pressure pad sensor according to the need of use (N.W indicates that both sides are off, Single represents only the pre-warning area near one end of the transmitter, Double represents the pre-warning area at both ends of the pressure pad, when the human gravity pressure is placed in the pre-warning area, the alarm main unit will immediately sound alarm).
- Leave bed testing:  
a. simulate the human body lying or sitting on the pressure pad and then leave, the pressure pad sensor immediately sends a signal (Model 636 pressure pad sensor delay 3 seconds), activates the alarm to sound alarm or vibration while

**PRODUCT PARAMETERS**

Model	Cushion size	thickness	color
636	10*35inch+14*14inch	0.2inch	White+black
636S	10*35inch	0.2inch	White

**CORRECT DISPOSAL**  
This product is subject to the provision of European Directive 2012/19/EU. The symbol showing a wheeled bin crossed through indicates that the product requires separate refuse collection in the European Union. This applies to the product and all accessories marked with this symbol. Products marked as such may not be discarded with normal domestic waste, but must be taken to a collection point for recycling electrical and electronic devices.

**FCC INFORMATION**  
**CAUTION:** Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment! This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:  
1) This product may cause harmful interference.  
2) This product must accept any interference received, including interference that may cause undesired operation.

**WARNING:** Changes or modifications to this product not expressly approved by the party responsible for compliance could void the user's authority to operate the product.

**Note:** This product has been tested and found to comply with the limits for a Class B digital device pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This product generates, uses and can radiate radio frequency energy, and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that

**VEVOR**  
TOUGH TOOLS, HALF PRICE  
Technical Support and E-Warranty Certificate  
[www.vevor.com/support](http://www.vevor.com/support)

**Fall Mats**  
Model: 636, 636S

We continue to be committed to provide you tools with competitive price. "Save Half", "Half Price" or any other similar expressions used by us only represents an estimate of savings you might benefit from buying certain tools with us compared to the major top brands and does not necessarily mean to cover all categories of tools offered by us. You are kindly reminded to verify carefully when you are placing an order with us if you are actually Saving Half in comparison with the top major brands.

210mm

the sensor area of the sensor LED (Go on forever). Press the Reset key twice to release the current alarm, and the main unit enters the working standby state again. The wireless distance between the alarm main unit and the pressure pad sensor is about 40 to 50 meters.

**CAUTIONS**

- The surface material of the pressure pad is PVC. Please don't direct touch with human skin when using it. You can put the pressure pad under the bed sheet or thin cushion and test whether it can work normally.
- Alarm main unit. The red LED indication when the pressure pad sensor battery is low, and the LED light dims when the button battery is low. Please replace the battery in time to ensure normal use.
- If the product is not used for a long time, please take out the battery to avoid battery leakage and damage to the device.
- This product is an electronic safety auxiliary product, and users must confirm the necessary reliability according to the usage properties.
- If you have any questions during the use of the product, please contact with local distributor.

**Pre-warning testing (For 636):** The simulated human body lies on the pressure pad sideways or moves to the warning area at both ends of the pressure pad placed on the edge of the bed. In that case, the pressure pad sensor detects that the human body will falling from the bed is in danger and immediately sends a signal to activate the alarm main unit to prompt family members or caregivers to adjust the bed position of the caregiver to the safe area in the center of the bed.

**Emergency CALL test:** press the call button and activate the alarm host to prompt family members or caregivers to help.

**Alarm mode/volume setting:** Tap the MODE key repeatedly to set the mode and volume you want, high volume + vibration/Low volume + Vibration/Vibration mode only.

**Learn code pairing and deletion**

- The product has learn code paired when left factory, if you need to change the pairing SENSOR area or increase the pairing use (main unit can be paired with 4 sensors), long press the LEARN key for 5 seconds to enter the pairing mode, SENSOR 1 LED flashing, indicating that the area can be paired. Trigger the SENSOR to be paired (the pressure pad SENSOR can press the pressure pad with the palm and then loose the palm, and the CALL button can press the CALL key), and hear the Bi sound of the main unit. SENSOR 1 If the LED is long on, the pairing is successful. For example, if the SENSOR 1 LED is long on, it indicates that there is pairing in the area. Tap the LEARN key to select the SENSOR 1-4 pairing area and repeat the SENSOR 1 pairing operation for pairing. Press the RESET key or no operation for 30 seconds to exit the learning code pairing mode.
- Delete the SENSOR pairing. At the same time, Press the MODE+RESET key for 5 seconds to enter the deletion mode and the SENSOR 1 LED will light long on/press the LEARN key to recycle select the existing SENSOR 1-4 (for deletion), and hold down the LEARN key for 3 seconds to delete the pairing in the SENSOR 1 area. SENSOR 1 LED Blinking indicates that the SENSOR 1 is successfully deleted. SENSOR 2-4 can repeat SENSOR 1 operation. Press the RESET key or do not perform any operation for 30 seconds to exit the deletion mode.

**UK REP**  
YH CONSULTING LIMITED.  
CID YH Consulting Limited Office 147,  
Centurion House, London Road,  
Staines-upon-Thames, Surrey, TW18 4AX

**EC REP**  
E-CrossStu GmbH  
Mainzer Landstr.69,  
60329 Frankfurt am Main.

**VEVOR**  
TOUGH TOOLS, HALF PRICE  
Technical Support and E-Warranty Certificate  
[www.vevor.com/support](http://www.vevor.com/support)

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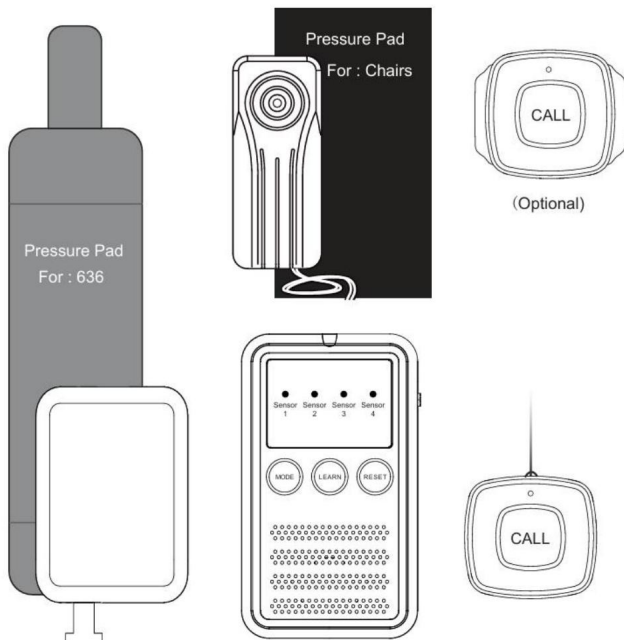
[www.vevor.com/support](http://www.vevor.com/support)

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**Importerad till AUS:** SIHAO PTY LTD. 1 ROKEVA STREETEASTWOOD NSW 2122
**Importerad till USA:** Sanven Technology Ltd. Suite 250, 9166 Anaheim Place, Rancho
C/O YH Consulting Limited Office 147,