VEVOR[®]

TOUGH TOOLS, HALF PRICE

Technical Support and E-Warranty Certificate
www.vevor.com/support

Red Light Therapy Wand

MODEL:TL09-B

We continue to be committed to provide you tools with competitive price. "Save Half", "Half Price" or any other similar expressions used by us only represents an estimate of savings you might benefit from buying certain tools with us compared to the major top brands and does not necessarily mean to cover all categories of tools offered by us. You are kindly reminded to verify carefully when you are placing an order with us if you are actually saving half in comparison with the top major brands.

Model: TL09-B

NEED HELP? CONTACT US!

Have product questions? Need technical support? Please feel free to contact us:

Technical Support and E-Warranty Certificate
www.vevor.com/support

This is the original instruction, please read all manual instructions carefully before operating. VEVOR reserves a clear interpretation of our user manual. The appearance of the product shall be subject to the product you received. Please forgive us that we won't inform you again if there are any technology or software updates on our product.

Congratulations on taking control of your pain !
Now you can use the LED Red Light Therapy to provide direct, non-addictive, drug free pain relief.

Efficient Cells Can Enhance Health

Uses a combo of red 630nm 660nm and Infrared 710nm 850nm 900nm LEDs which is chosen to give the ideal light intensity for biological effects anywhere on the body.

1. Included With This Package

1*Red Light Therapy Wand

1*Carry Bag

2*Rechargeable Battery18650

1*Battery Charger

- 1*USB Cable
- 1*Protective glasses
- 1*User Manual
- 1*Decorative ring

2. Benefits

Temporary relief of minor muscle and joints pain, arthritis and muscle spasm, relieving stiffness, promoting relaxation of muscle tissue.

3. Directions

1)Commit to a 3-week Treatment Plan.While many people achieve pain relief from Red Light Therapy almost immediately,others see measurable pain relief incrementally over time.Depending upon injury type severity and different pain pathways involved,people vary in therapy response time.

This is why we recommend committing to a 3-week Treatment Plan. If you see results right away,great,continue treating as needed. However, if you do not feel relief as quickly, know that this is also normal, so try not to get frustrated Simply treat your pain onsistently with 3 quick treatments per day with the, For up to 3 minutes each time per treatment spot. After 3 weeks of consistent red light therapy, most people experience a diminishing of their pain,

2) Do not suddenly stop your other treatments upon beginning red light therapy. Although you may not be satisfied with current

treatments, they may be having a partial effect on the level of your pain. While continuing with your current treatments, start your red light therapy, giving your red device time to provide optimal pain relief. After you have achieved maximum pain relief from your red device, you may then consider slowly weaning yourself off of other medications or treatments. If you plan to reduce or stop a prescription medication, first consult with your physician.

For many people, it works so well that they are able to come off pain pills completely. Others are able to reduce the amount of pain medications they are on, reducing the risk of drug side effects.

Hold the unit over the painful area at a distance anywhere from touching the skin to 1 inch away from the skin, depending on the depth of treatment desired. The Red light penetrates up to 1 inch into your tissues or joints.

Pain Level	Light Pain	Average pain	Deep pain
DOSAGE	1 DOSE Per treatment spot (1 minute treatment) Once a day	2 DOSES Per treatment spot 2 minutes treatment Twice or 3x a day	3 DOSES Per treatment spot 13 minutes treatment 3 times a day

Do not keep the device on for more than 5 minutes at a time. If you want to treat several parts of your body consecutively, turn off the device for at least 1 minute, after every 5 minutes on. This will keep the unit from overheating.

It has two using ways, one you can use it on your mouth, nose, ear.

Another way you can remove the Adapter of device, and you can use it on your wrist, finger, knee.

		
For mouth treatment	For ear treatment	For nose treatment
		
For wrist treatment	For finger treatment	For knee treatment

How to install batteries?

Insert the Battery provided on the battery compartment. -Positive Pole against the TOP.

-Negative Pole against the BOTTOM.

*Please only use 18650 batteries.

*18650 batteries can be charged for 2-3 hours.

TOP 

BOTTOM 

4. Product Specification

Product Size	φ29.5x H 216mm	Net Weight	0.4kg
Materia	Aluminium alloy	Wavelength	630:660nm 710:850:900 nm
Irradiance	Surface:>180mw/cm2	Deiault time	5 minutes
Powe	5.4W	Battery Capacity	2500mAh
Warranty	1 year	Lifespan	50,000 Hours

Function Instruction

1. Wavelength: 630nm, 660nm, 710nm, 850nm, 900nm

2. Timer: 5minutes default

3. Three modes:

(a) Mode1 (All LEDs working):
630nm, 660nm, 710nm, 850nm, 900nm

(b) Mode2 (Only red working):
630nm, 660nm;

(c) Mode3 (Only Near infrared working):
710nm, 850nm, 900nm;

4. Support pulsing function.

5. Button operation:

(a) Press first time: Mode1;

(b) Press second time: Mode2;


(c) Press third time: Mode3;

(d) Press fourth time: Mode1+10Hz Pulsing;

(e) Press fifth time: Mode2+10Hz Pulsing;

(f) Press sixth time: Mode3+10Hz Pulsing;

(g) Press seventh time: off.

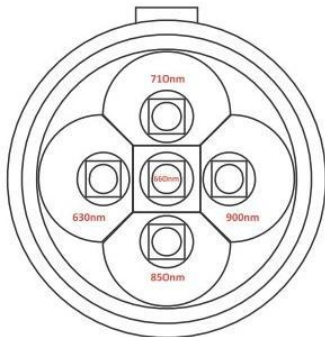Research points to a brief "quench period" that takes place in the cells when NIR light is pulsed. Although it lasts only a few milliseconds, this short pause makes a big difference, allowing cells that are under oxidative stress to respond even better to

light. Then it is much better for cells stimulation and recovery with pulsing therapy.

CORRECT DISPOSAL

This product is subject to the provision of European Directive 2012/19/EU. The symbol showing a wheeled bin crossed through indicates that the product requires separate refuse collection in the European Union. This applies to the product and all accessories marked with this symbol. Products marked as such may not be discarded with normal domestic waste, but must be taken to a collection point for recycling electrical and electronic devices.

BATTERY DISPOSAL

To minimize hazards to health and the environment at the end of this product's life, laws dealing with Waste Electrical and Electronic Equipment (WEEE) and The Waste Battery Directive require you to dispose of this product at a suitable collection facility where it will be sent in order to remove the batteries and for appropriate recycling. Please contact your local authorities for more details on recycling and safe disposal of these in your area.

FCC Information:

CAUTION: Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment!

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

1) This product may cause harmful interference.

2) This product must accept any interference received, including interference that may cause undesired operation.

WARNING: Changes or modifications to this product not expressly approved by the party responsible for compliance could void the user's authority to operate the product.

Note: This product has been tested and found to comply with the limits for a Class B digital device pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This product generates, uses and can radiate radio frequency energy, and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this product does cause harmful interference to radio or television reception, which can be determined by turning the product off and on, the user is encouraged to try to correct the interference by one or more of the following measures.

- Reorient or relocate the receiving antenna.
- Increase the distance between the product and receiver.
- Connect the product to an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for assistance.

5. Caution

- *Don't use the light around water
- *Don'tstare directly at the light use protective glasses while using the device
- *Don'ttreat near the eyes

6. Equipment Maintenance

Cleaning And Disinfecting:

Use only a soft dry cloth to clean the exterior. Never use an abrasive cleanser or submerge in water or any other liquid. This device cannot be sterilized by any liquid or autoclave method.

Storage: Store the unit in a cool, dry place, away from direct light and extreme heat.

* Warning: Do not allow moisture to come in contact with the case or light of the unit. Do not immerse the unit in water.

* Warning: If unit will not be used for an extended period of time (1 month+) remove the batteries from the unit to avoid their Corrosion.

VEVOR[®]

TOUGH TOOLS, HALF PRICE

Technical Support and E-Warranty Certificate

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