

# **VEVOR<sup>®</sup>**

**TOUGH TOOLS, HALF PRICE**

Technical Support and E-Warranty Certificate [www.vevor.com/support](http://www.vevor.com/support)

## **RED LIGHT THERAPY BELT**

**MODEL: BRL-1120**

We continue to be committed to provide you tools with competitive price.


"Save Half", "Half Price" or any other similar expressions used by us only represents an estimate of savings you might benefit from buying certain tools with us compared to the major top brands and does not necessarily mean to cover all categories of tools offered by us. You are kindly reminded to verify carefully when you are placing an order with us if you are actually saving half in comparison with the top major brands.

# VEVOR®

TOUGH TOOLS, HALF PRICE

## RED LIGHT THERAPY BELT

MODEL: BRL-1120

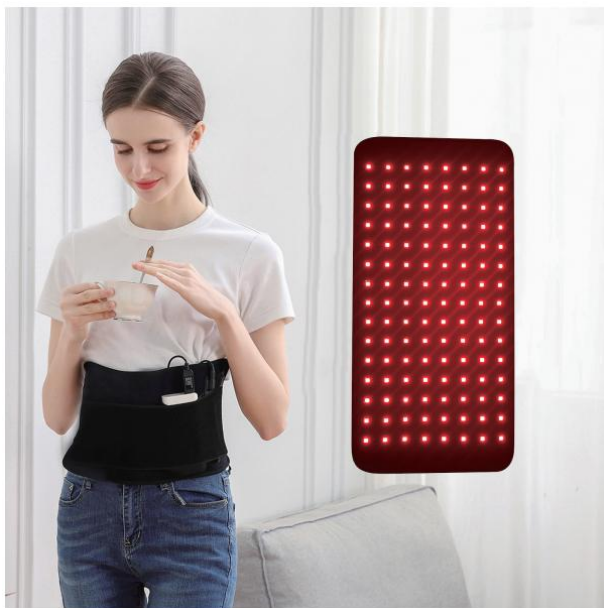

### NEED HELP? CONTACT US!

Have product questions? Need technical support? Please feel free to contact us:

**Technical Support and E-Warranty Certificate**  
**[www.vevor.com/support](http://www.vevor.com/support)**

This is the original instruction, please read all manual instructions carefully before operating. VEVOR reserves a clear interpretation of our user manual. The appearance of the product shall be subject to the product you received. Please forgive us that we won't inform you again if there are any technology or software updates on our product.

# Instructional Manual

## Red-Light-Therapy for Human & Animals!

Thank you for purchasing our light, we hope you can reap benefits of red light therapy at home.

### Say hello to your red light!

Before starting your first session, please make sure to read this manual as it may explain important operating procedures.

**Notes:** Each led diode inside with three chips at 1 pc 660nm and 2 pcs 850nm is invisible to the human made eyes.

## Specifications

Rating Power	60W	Wavelength	1*660nm 2*850nm
Actual Power	15W	Irradiance	Max: 40m W/c m <sup>2</sup>
Product Size	16.1in*7.9in/41cm*20cm	Timer	20min
LED QTY	120*0.5w	Material	Fabric
Net Weight	200g	Lifespan	Over 50,000 Hours

The red light device uses combo red 660nm and Infrared 850nm leds which is chosen to give the ideal wavelengths found in the scientific research of the last three decades.

660nm of red light

This light can be seen by the human eyes and targets skin cells.

850nm of near infrared light

This light cannot be seen by the human eyes and targets deep tissue.

This package includes

- 1\* Red light therapy pad
- 1\* Adaptor
- 1\* Stretchy straps

- 1\* Power cord
- 1\* USB power cord
- 1\* User manual

Way 2 : Power bank supply power



**Benefits backed by science:**

Research shows that the wavelengths in the spectrum of 660-900nm show optimalbiotical response.

- Muscle and Joints Pain
- Speed wound healing
- Enhanced blood circulatorn
- Weight loss
- Anti-inflammatory
- Decreases pain

**Even though this technology don' t have side effects , we suggest taking the following precautions:**

This device is not waterproof

This device is not suitable for children and pregnant woman

Do not use this device exceed 30 minutes each time

Do not stare at the light directly

It's NOT recommend to use it when sleeping


Please storage it away from children

## Application

Here show you some examples how to use the red light device . You and your pets can enjoy the red light at different ways for like thighs, calves, knees.

Near Infrared-Infrared A(IR A)

Usually use 660nm~850nm

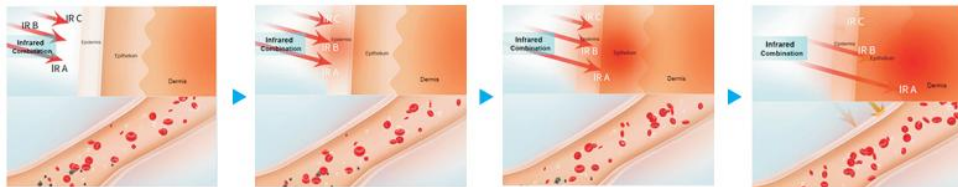

Invisible near infrared energy are absorbed by muscle , subcutaneous tissue, which results in THERMAL EFFECT, your body responds to this radiation heat with increased blood flow,allowing the parts of the body to receive the oxygen and nutrients they need in order to function more effectively,regeneration is stimulated.

## Measured value



### Note:

**Indoor use only :** Light is not waterproof . can not used in water or in wet places. Especially can not be used with essential oils.

**Repairing :** If your light need after-sale service , please contact us directly.

**Installation:**To avoid the danger of electric shock , please turn off the power switch before installation.



## SAFETY CAUTION:

1. Do not make direct view or observation to the laser beams via any optical instruments.
2. Keep out of reach of children to prevent injuries to their eyes.
3. It's forbidden for pacemaker users since the device may interfere with the operation of pacemaker.
4. During the treatment period, the diet structure should be improved and proper exercises are necessary to enhance the metabolism.
5. When using the device, please keep light emitting ports clean. Otherwise, the light losses will increase and then affect effects.
6. The symptoms like thirst, dry eyes, dizziness or rise of pulse rate after use are normal phenomena.

## CORRECT DISPOSAL



This product is subject to the provision of European Directive 2012/19/EU. The symbol showing a wheellie bin crossed through indicates that the product requires separate refuse collection in the European Union. This applies to the product and all accessories marked with this symbol. Products marked as such may not be discarded with normal domestic waste, but must be taken to a collection point for recycling electrical and electronic devices.

## FCC Information:

**CAUTION:** Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment!

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- 1) This product may cause harmful interference.
- 2) This product must accept any interference received, including interference that may cause undesired operation.

**WARNING:** Changes or modifications to this product not expressly approved by the party responsible for compliance could void the user's authority to operate the product.

**Note:** This product has been tested and found to comply with the limits for a Class B digital device pursuant to Part 15 of the FCC Rules, These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This product generates, uses and can radiate radio frequency energy, and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this product does cause harmful interference to radio or television reception, which can be determined by turning the product off and on, the user is encouraged to try to correct the interference by one or more of the following measures.

- Reorient or relocate the receiving antenna.
- Increase the distance between the product and receiver.
- Connect the product to an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for assistance.

**VEVOR<sup>®</sup>**

**TOUGH TOOLS, HALF PRICE**

**Technical Support and E-Warranty Certificate**

**[www.vevor.com/support](http://www.vevor.com/support)**