VEVOR

Affordable. Reliable. Home Improvement.

LAB DISTILLATION APPARATUS


USER MANUAL

MODEL:YHEO-7

VEVOR
Affordable. Reliable. Home Improvement.

LAB DISTILLATION APPARATUS

MODEL:YHEO-7

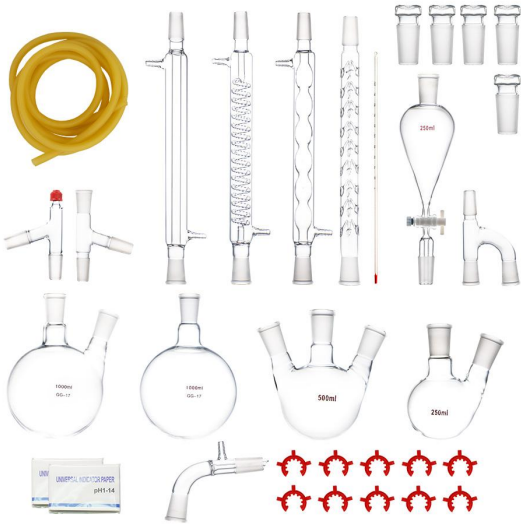This is the original instruction, please read all manual instructions carefully before operating. VEVOR reserves a clear interpretation of our user manual. The appearance of the product shall be subject to the product you received. Please forgive us that we won't inform you again if there are any technology or software updates on our product.

SAFETY INSTRUCTIONS

WARNING:

Read this material before using this product. Failure to do so can result in serious injury.

Assembly precautions

1. Assemble only according to these instructions. Improper assembly can create hazards.
2. Wear safety goggles and duty work gloves during assembly.
3. Keep the assembly area clean and well-lit.
4. Keep bystanders out of the area during assembly.
5. Do not assemble when tired or when under the influence of alcohol, drugs or medication.
6. The product capabilities apply to properly and completely assembled products only.
7. Assemble on a flat, level, hard and smooth surface capable of safely supporting the Lab Distillation Apparatus.
8. For additional information regarding the parts listed in the following pages, please refer to the Assembly Diagram of this manual. Unwrap and separate all parts in a clean work area.

Use precautions


1. DO NOT SIT OR STAND ON THIS ITEM.
2. This product is not a toy. Do not allow children to play with or near this item.
3. Do not exceed weight capacities.
4. Use only on a flat, level, hard and smooth surface capable of safely supporting a fully loaded Lab Distillation Apparatus.
5. Use as intended only.
6. Inspect before every use; do not use if parts are loose or damaged.

SAVE THIS MANUAL

INSTRUCTIONS

Thank you very much for choosing this model.Glass fittings should always be handled with care. Please read all of the instructions before using it.The information will help you achieve the best possible results.

Parts name	Number of parts	Description
flask	1	1000ml,24/40 ground glass joint
Two-necked flask	1	250ml,24/40 ground glass joint
Two-necked flask	1	1000ml,24/40 ground glass joint
Three-necked flask	1	500ml,24/40 ground glass joint
Straight condenser	1	/
Coil condenser	1	/
Ball condenser	1	/
Vigreaux column	1	/
Vacuum adapter	1	/
Distilling head	1	/
Y-tube	1	/
Thermometer adapter	1	/
Separatory funnel	1	250ml
Glass stopper	5	/
Thermometer	1	200℃
Plastic clip	10	/
Latex tubing	1	2m
PH test paper	2	/

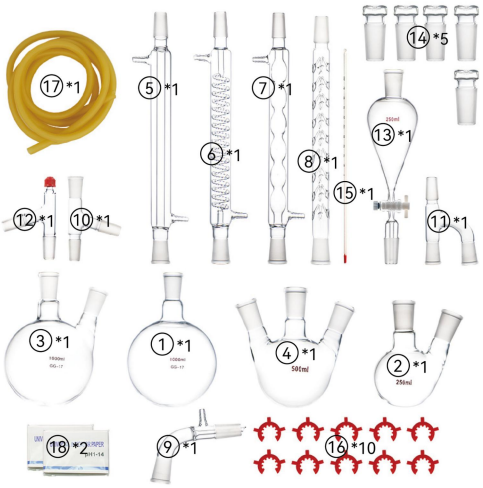No.	Parts Name
1	Flask
2	Two-necked flask
3	Two-necked flask
4	Three-necked flask
5	Straight condenser
6	Coil condenser
7	Ball condenser
8	Vigreaux column
9	Vacuum adapter
10	Distilling head
11	Y-tube
12	Thermometer adapter
13	Separatory funnel
14	Glass stopper
15	Thermometer
16	Plastic clip
17	Latex tubing
18	PH test paper

DIVERSE APPLICATION

Distillation apparatus kit is suitable for chemistry, medical treatment industries, also great for teaching activity, DIY essential oil.

PRODUCT PARAMETER

Model	YHEO-7
Number of Kits	32
Glass Material	Borosilicate Glass

Manufacturer: Shanghaimuxinmuyeyouxiangongsi**Address:** Shuangchenglu 803nong11hao1602A-1609shi, baoshanqu, shanghai 200000 CN.
Imported to AUS: SIHAO PTY LTD. 1 ROKEVA STREETEASTWOOD NSW 2122 Australia
Imported to USA: Sanven Technology Ltd. Suite 250, 9166 Anaheim Place, Rancho Cucamonga, CA 91730

UK	REP

YH CONSULTING LIMITED. C/O YH Consulting Limited Office 147, Centurion House, London Road, Staines-upon-Thames, Surrey, TW18 4AX

EC	REP

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Importerad till Australien: SIHAO PTY LTD. 1 ROKEVA STREETEASTWOOD NSW 2122 Australien
Importerad till USA: Sanven Technology Ltd. Suite 250, 9166 Anaheim Place, Rancho Cucamonga, CA 91730