VEVOR


Upgrade · The Home Creator Way

Tub Transfer Benches

MODEL: SC6045C

Tub Transfer Benches

MODEL:SC6045C

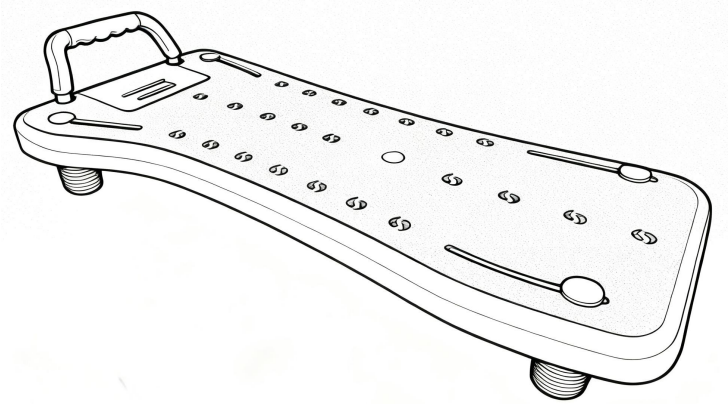This is the original instruction, please read all manual instructions carefully before operating. VEVOR reserves a clear interpretation of our user manual. The appearance of the product shall be subject to the product you received. Please forgive us that we won't inform you again if there are any technology or software updates on our product.

WARNING: Read and understand this entire manual before operating or servicing this product. Failure to follow these warnings and instructions can cause personal injury or damage to valuable property.

- Avoid children using the product. And this product is not a toy. Do not allow children to play.

ASSEMBLY PRECAUTIONS

1. Assemble only according to these instructions. Improper assembly can create hazards.
2. Wear ANSI-approved safety goggles and heavy-duty work gloves during assembly and use.

3. Keep the assembly area clean and well-lit.
4. Keep bystanders out of the area during assembly.
5. Do not assemble when tired or when under the influence of alcohol, drugs or medication.
6. Product capabilities apply to properly and completely assembled product only.
7. The product complies with relevant standards and regulations.Installation of this device must comply with all applicable regulations,and it should only be operated in well-ventilated areas.
8. This product is not intended for use by persons (including children) with reduced physical,sensory or mental capabilities or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
9. Please properly keep product.Check whether the product is abnormal before use. If it is abnormal,please stop using it.
10. Please use the product according to the intended use, warn of potential injury from misuse.
11. For additional information regarding the parts listed in the following pages, please refer to the Assembly Diagram of this manual. Unwrap and separate all parts in a clean work area. Please keep small spare parts out of children's reach.

SAVE THESE INSTRUCTIONS

WARNING:



Read the manual completely to avoid harm, product damage, or warranty avoidance. Seek professional help if you encounter discomfort, pain or worsening symptoms. If necessary, use supervision in order to assemble and use this product correctly and safely. Do not exceed the weight limit.
Regularly check for secure assembly and contact support if in doubt. Ensure the product’s legs are even and locked before use, and discontinue usage if there are signs of wear, damage, or missing parts. The product is to be used on flat, stable surfaces. Do not disassemble if not necessary, as it might compromise the longevity of the product. Do not modify the product unless specified in the manual.

IMPORTANT

1. Don’t exceed the maximum weight capacity.
2. Only for use on bathtubs that are stable and have flat surfaces on both sides.
3. Regularly inspect for signs of wear, damage, or loose parts.
4. Clean regularly using a mild detergent and warm water, and rinse off any soap residue. After cleaning, dry thoroughly.
5. Before use, confirm with your healthcare provider that the product meets your specific needs. The following guidelines are not a substitute for medical advice.

PART LIST

When unpacking this product, check to make sure the parts listed below are included and carefully inspect them for any damage that may have occurred during transit. If any part is missing or damaged, do not attempt to assemble or use the product.

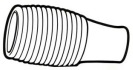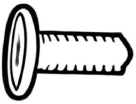NO	NAME	PICTURE	QTY
1	Handrail		1
2	Fixed leg		4
3	Screw A		4
4	Washer		4
5	Screw B		2
6	Manual	/	1


PRODUCT PAREMETERS

Model	SC6045C
Main Materials	PP
The applicable inner width of the bathtub	41-62cm
Seat Rated Load	150kg/330lbs


INSTALLATION AND OPERATION INSTRUCTIONS

1. Installation Procedures


As shown in the figure, the various parts of the bath chair are combined. The armrests are fixed in the two installation holes at the end through screw B, and 4 fixed legs are respectively fixed in the four slots of the seat board through screw A.

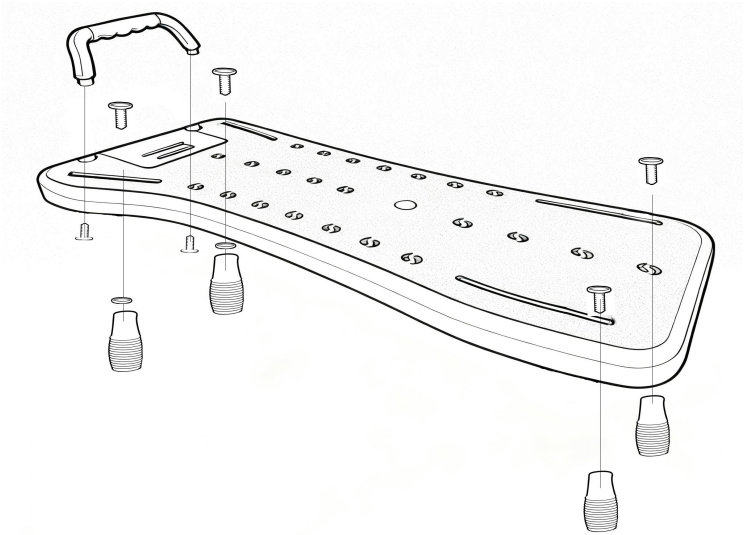
2. Adjustment Procedures

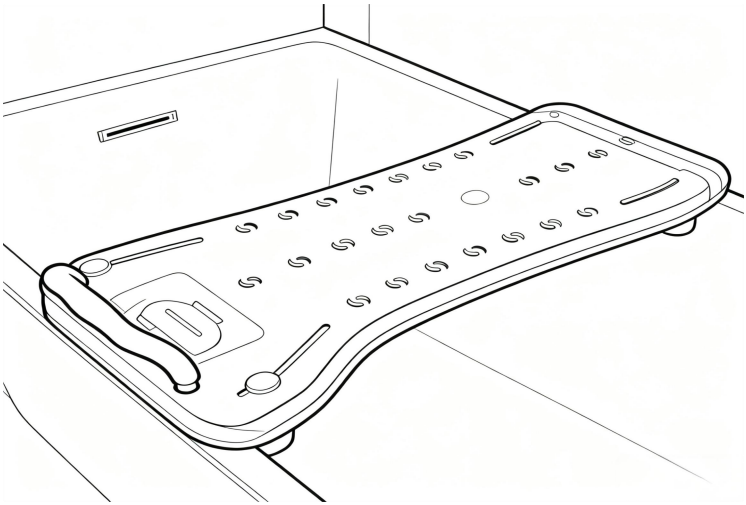Press the fixed legs tightly against the bathtub wall, and then tighten them to fix the seat plate in the width direction of the bathtub.

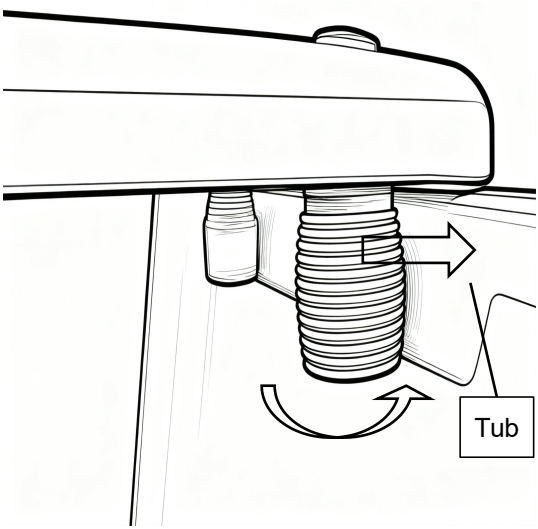The 4 fixed legs can move left and right to fit the applicable bathtub width range, which is 41-62cm.

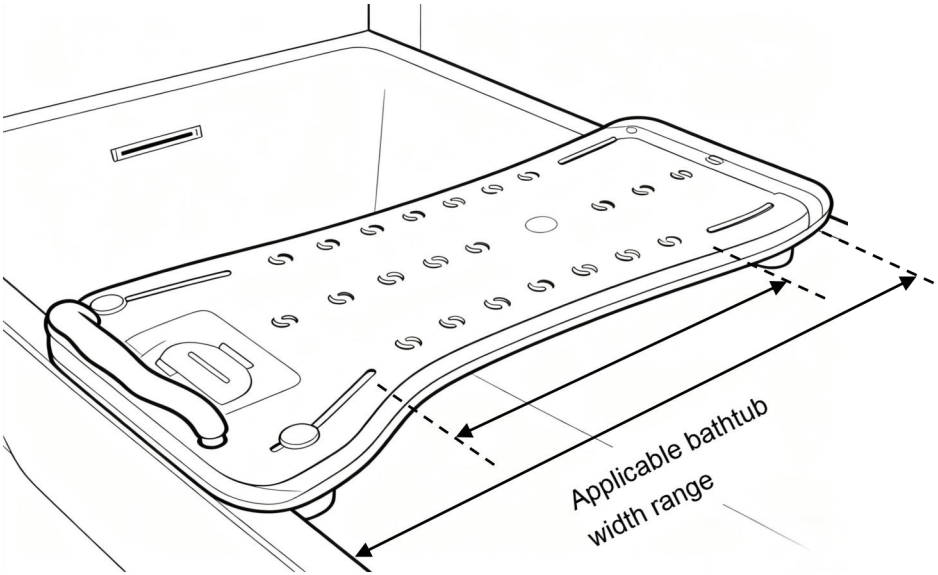Manufacturer: FOSHAN SUNCARE MEDICAL PRODUCTS CO..LTD
Address: 64-65 Building, zhongnan Industry Garden,No.88Haitian Road, Gaoming District, Foshan City, Guangdong Province, P.R.China
amy@suncare-medical.com
Imported to AUS: SIHAO PTY LTD. 1 ROKEVA STREETEASTWOOD NSW 2122 Australia
Imported to USA: Sanven Technology Ltd. Suite 250, 9166 Anaheim Place, Rancho Cucamonga, CA 91730

EC REP

INDEPSALE TECHNOLOGY LIMITED
Office 403, Nesta Business Centre, 4-5 Burton Hall Rd,Sandyford, Dublin 18, D18 A094, Ireland
Email: support@indepsale.ie Tel:+353874373999

UK REP

VEVOR SMART LIMITED
2 Goldstone Close,Brighton and Hove, Hove, BN3 7PD, UnitedKingdom
Email: supports@vevor.co.uk Tel:+447795204660

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Tillverkare: FOSHAN SUNCARE MEDICAL PRODUCTS CO.,LTD
Importerad till Australien: SIHAO PTY LTD. 1 ROKEVA STREETEASTWOOD NSW 2122 Australien
Importerad till USA: Sanven Technology Ltd. Suite 250, 9166 Anaheim Place, Rancho Cucamonga, CA 91730