VEVOR

Affordable. Reliable. Home Improvement.

Raised Toilet Seat

Model: SC7060H

VEVOR
Affordable. Reliable. Home Improvement.

Raised Toilet Seat

Model: SC7060H


This is the original instruction, please read all manual instructions carefully before operating. VEVOR reserves a clear interpretation of our user manual. The appearance of the product shall be subject to the product you received. Please forgive us that we won't inform you again if there are any technology or software updates on our product.


Thank you very much for choosing this product. Please always handle it with care. Please read all of the instructions before using it. The information will help you achieve the best possible results.

SPECIFICATION

Model	Raised height	load
SC7060H	50/100/150mm	300lbs/136kg

- Adapt to toilet： Toilet interior width ≥ 240mm, toilet exterior width ≤ 350mm, suitable depth of the toilet: from 430mm to 470mm.
- Please check the size of your toilet before choosing the raised toilet seat.

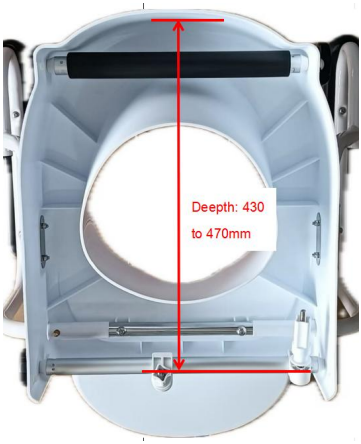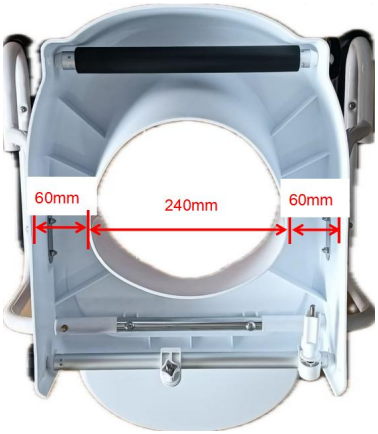PACKAGE LIST

No.	Name	Picture	**SC7060H**
1	Raised Toilet Seat with cover		1
2	Armrest		2

3	Manual		1

CAUTIOUS

1. Always be sure that the seat is correctly and securely locked in place before using.
2. Be sure that the locking of the seat is tight at all times, thus securing the seat to the toilet bowl. Physically challenged users should be supervised or assisted whenever using this product.
3. The raised toilet seat MUST NOT be used when overloaded. User's weight should be centered over toilet seat, not too forward or too far to one side, as unit may tip.
4. This product is designed for use with the most common types of commodes. It will not fit all commode types; be sure that the unit is correctly installed, fits properly and is stable before use.
5. For models with armrests: The arms are for the purpose of providing assistance when getting on and off the seat. DO NOT try to use the arms to support full body weight during transfer, as this may result in failure of the unit and risk of injury.
6. If you have a style of commode that this product will not fit properly, or cannot be tightened securely as noted above, it is recommended that you discontinue usage.
7. Regularly check for secure assembly and contact support if in doubt. Ensure the product's legs are even and locked before use, and discontinue usage if there are signs of wear, damage, or missing parts. Do not disassemble if not necessary, as it might compromise the longevity of the product. Do not modify the product unless specified in the manual.
8. Exercise caution during operation and assembly. Avoid extreme temperatures or open flames. If an allergic reaction occurs, cease use immediately and consult a healthcare professional.
9. Before use, confirm with your healthcare provider that the product meets your specific needs. The following guidelines are not a substitute for

medical advice.

INSTRUCTIONS


1. Remove your existing toilet seat from the bowl. Seats are generally attached by a pair of bolts/nuts located under the rear rim of the toilet.
Unscrew and remove the long plastic nuts from the rear of the toilet seat raiser.


2. Position the unit onto your toilet bowl, inserting the rear mounting bolts through the existing holes on your toilet.
Ensure the unit is centred over the toilet bowl.

3. The position of the mounting bolts can be adjusted to match different toilets. To adjust, slide the mounting bolts along the rear bar to achieve the correct distance.

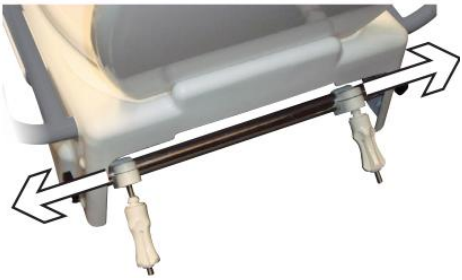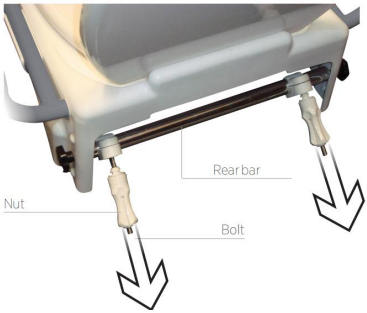Mounting nuts must be loose for mounting bolts to slide freely along the bar.

Reattach both long plastic nuts to the mounting bolts and tighten firmly.
4. Check that the toilet seat raiser is firmly clamped against the toilet bowl before use.

Adjusting seat height

The toilet seat raiser can be adjusted to 3different height settings, adding 50,100 or 150 mm to the height of the toilet bowl Height can be adjusted either before or after installation.

When adjusting after the toilet seat raiser is installed, be careful not to drop components into the toilet.

Armrests assembly

The toilet seat raiser can be used with or without armrests.One or both armrests can be used to suit individual users and/or bathroom set up. A brass locking button is used to secure each armrest.

To insert the armrests, press in the locking button located under the frame

of the armrest and slide the armrest into the arm hole at the side of the toilet seat. Check that the armrest is stable and secure before use.

MAINTENANCE

1. REGULAR CLEANING IS RECOMMENDED, USE A NON-ABRASIVE CLEANER OR DETERGENT. RINSE WITH WARM WATER AND DRY COMPLETELY. ABRASIVE CLOTHS MAY DAMAGE THE SURFACE.
2. PLEASE REGULARLY INSPECT YOUR RAISED TOILET SEAT WITH ARMS AND ITS PARTS FOR ANY CORROSION, BREAKAGE, COMPRESSION OR WEAR.
3. Regular maintenance

Regularly check that the mounting screws are tight and that the toilet seat raiser is firmly clamped to the toilet bowl.

Care & cleaning

To clean the unit after each use, wipe with a mild soap and warm water solution and rinse well. Dry with a clean cloth.

For heavier soiling, wipe clean with a diluted bleach solution and rinse well before drying.

Do not use water temperatures in excess of 80°c.

Do not use high pressure cleaners, steam cleaners or harsh, abrasive cleaning agents or utensils (such as a wire brush) as these can lead to deterioration of the unit.

Manufacturer: FOSHAN SUNCARE MEDICAL PRODUCTS CO., LTD.
Address: 2-3 Floor of 65 Building, Zhongnan Industry Garden, No.88 Haitian Road, Gaoming District, Foshan City, Guangdong, 528500, CHINA.
Imported to AUS: SIHAO PTY LTD. 1 ROKEVA STREETEASTWOOD NSW 2122 Australia
Imported to USA: Sanven Technology Ltd. Suite 250, 9166 Anaheim Place, Rancho Cucamonga, CA 91730

UK	REP

YH CONSULTING LIMITED. C/O YH Consulting Limited Office 147, Centurion House, London Road, Staines-upon-Thames, Surrey, TW18 4AX

EC	REP

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Tillverkare: FOSHAN SUNCARE MEDICAL PRODUCTS CO ., LTD .
Importerad till Australien: SIHAO PTY LTD. 1 ROKEVA STREETEASTWOOD
Importerad till USA: Sanven Technology Ltd. Suite 250, 9166 Anaheim Plats, Rancho Cucamonga, Kalifornien 91730
E-CrossStu GmbH